This Page Is Inserted by IFW Operations and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

As rescanning documents will not correct images, please do not report the images to the Image Problem Mailbox.

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification ⁶ : C12Q 1/68, G06F 15/00	A1	11) International Publication Number: WO 95/2068: 43) International Publication Date: 3 August 1995 (03.08.95
(21) International Application Number: PCT/US (22) International Filing Date: 27 January 1995 ((30) Priority Data: 08/187,530 27 January 1994 (27.01.94) 08/282,955 29 July 1994 (29.07.94) (71) Applicant: INCYTE PHARMACEUTICALS, INC. 3330 Hillview Avenue, Palo Alto, CA 94304 (US) (72) Inventors: SEILHAMER, Jeffrey, J.; 12555 La Cr Altos Hills, CA 94022 (US). SCOTT, Randal, W Sun-Mor, Mountain View, CA 94040 (US). (74) Agents: CAGE, Kenneth, L. et al.; Willian Brinks Hof & Lione, 2000 K Street, N.W., Suite 200, Washin 20006-1809 (US).	27.01.9 [US/US] [US/US] [US/US] [US/US] [US/US] [US/US]	(81) Designated States: AM, AU, BB, BG, BR, BY, CA, CN, CZ EE, FI, GE, HU, IP, KG, KP, KR, KZ, LK, LR, LT, LV, MD, MG, MN, MX, NO, NZ, PL, RO, RU, SI, SK, TJ, TI UA, UZ, VN, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI pater (BF, BJ, CF, CG, CI, CM, GA, GN, ML, MR, NE, SN, TT TG), ARIPO patent (KE, MW, SD, SZ). Published With international search report.

(54) Title: COMPARATIVE GENE TRANSCRIPT ANALYSIS

(57) Abstract

A method and system for quantifying the relative abundance of gene transcripts in a biological specimen. One embodiment of the method generates high-throughput sequence-specific analysis of multiple RNAs or their corresponding cDNAs (gene transcript imaging analysis). Another embodiment of the method produces a gene transcript imaging analysis by the use of high-throughput cDNA sequence analysis. In addition, the gene transcript imaging can be used to detect or diagnose a particular biological state, disease, or condition which is correlated to the relative abundance of gene transcripts in a given cell or population of cells. The invention provides a method for comparing the gene transcript image analysis from two or more different biological specimens in order to distinguish between the two specimens and identify one or more genes which are differentially expressed between the two specimens.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
ΑŪ	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	Æ	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Келуа	RO	Romania
CA	Canada	KG	Kyrgystan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic	SD	Sudan
CG	Congo		of Korea	SE	Sweden
CH	Switzerland	KR	Republic of Korea	SI	Slovenia
CI	Côte d'Ivoire	KZ	Kazakhstan	SK	Slovakia
CM	Cameroon	u	Liechtenstein	SN	Senegal
CN	China	LK	Sri Lanka	TD	Chad
CS	Czechoslovakia	LU	Luxembourg	TG	Togo
cz	Czech Republic	LV	Latvia	TJ	Tajikistan
DE	Germany	MC	Monaco	TT	Trinidad and Tobago
DK	Denmark	MD	Republic of Moldova	UA	Ukraine
ES	Spain	MG	Madagascar	US	United States of America
FI	Finland	ML	Mali	UZ	Uzbekistan
FR	France	MN	Mongolia	VN	Viet Nam
GA ·	Gabon		•		

COMPARATIVE GENE TRANSCRIPT ANALYSIS

1. FIELD OF INVENTION

The present invention is in the field of molecular biology and computer science; more particularly, the

5 present invention describes methods of analyzing gene transcripts and diagnosing the genetic expression of cells and tissue.

2. BACKGROUND OF THE INVENTION

Until very recently, the history of molecular biology

10 has been written one gene at a time. Scientists have
observed the cell's physical changes, isolated mixtures
from the cell or its milieu, purified proteins, sequenced
proteins and therefrom constructed probes to look for the
corresponding gene.

Projects to sequence the billions of bases in the human genome. These projects typically begin with dividing the genome into large portions of chromosomes and then determining the sequences of these pieces, which are then analyzed for identity with known proteins or portions thereof, known as motifs. Unfortunately, the majority of genomic DNA does not encode proteins and though it is postulated to have some effect on the cell's ability to make protein, its relevance to medical applications is not understood at this time.

A third methodology involves sequencing only the transcripts encoding the cellular machinery actively involved in making protein, namely the mRNA. The advantage is that the cell has already edited out all the non-coding DNA, and it is relatively easy to identify the protein-coding portion of the RNA. The utility of this approach was not immediately obvious to genomic researchers. In fact, when cDNA sequencing was initially proposed, the method was roundly denounced by those committed to genomic sequencing. For example, the head of the U.S. Human Genome project discounted CDNA sequencing as not valuable and refused to approve funding of projects.

In this disclosure, we teach methods for analyzing DNA, including cDNA libraries. Based on our analyses and

research, we see each individual gene product as a "pixel" of information, which relates to the expression of that, and only that, gene. We teach herein, methods whereby the individual "pixels" of gene expression information can be combined into a single gene transcript "image," in which each of the individual genes can be visualized simultaneously and allowing relationships between the gene pixels to be easily visualized and understood.

We further teach a new method which we call electronic subtraction. Electronic subtraction will enable the gene researcher to turn a single image into a moving picture, one which describes the temporality or dynamics of gene expression, at the level of a cell or a whole tissue. It is that sense of "motion" of cellular machinery on the scale of a cell or organ which constitutes the new invention herein. This constitutes a new view into the process of living cell physiology and one which holds great promise to unveil and discover new therapeutic and diagnostic approaches in medicine.

We teach another method which we call "electronic northern," which tracks the expression of a single gene across many types of cells and tissues.

Nucleic acids (DNA and RNA) carry within their sequence the hereditary information and are therefore the prime molecules of life. Nucleic acids are found in all living organisms including bacteria, fungi, viruses, plants and animals. It is of interest to determine the relative abundance of different discrete nucleic acids in different cells, tissues and organisms over time under various conditions, treatments and regimes.

All dividing cells in the human body contain the same set of 23 pairs of chromosomes. It is estimated that these autosomal and sex chromosomes encode approximately 100,000 genes. The differences among different types of cells are believed to reflect the differential expression of the 100,000 or so genes. Fundamental questions of biology could be answered by understanding which genes are transcribed and knowing the relative abundance of transcripts in different cells.

Previously, the art has only provided for the analysis of a few known genes at a time by standard molecular bi logy techniques such as PCR, northern blot analysis, or other types of DNA probe analysis such as in situ

5 hybridization. Each of these methods allows one to analyze the transcription of only known genes and/or small numbers of genes at a time. Nucl. Acids Res. 19, 7097-7104 (1991); Nucl. Acids Res. 18, 4833-42 (1990); Nucl. Acids Res. 18, 2789-92 (1989); European J. Neuroscience 2, 1063-1073

10 (1990); Analytical Biochem. 187, 364-73 (1990); Genet. Annals Techn. Appl. 7, 64-70 (1990); GATA 8(4), 129-33 (1991); Proc. Natl. Acad. Sci. USA 85, 1696-1700 (1988); Nucl. Acids Res. 19, 1954 (1991); Proc. Natl. Acad. Sci. USA 88, 1943-47 (1991); Nucl. Acids Res. 19, 6123-27

15 (1991); Proc. Natl. Acad. Sci. USA 85, 5738-42 (1988); Nucl. Acids Res. 16, 10937 (1988).

Nucl. Acids Res. <u>16</u>, 10937 (1988).

Studies of the number and types of genes whose transcription is induced or otherwise regulated during cell processes such as activation, differentiation, aging, viral transformation, morphogenesis, and mitosis have been

pursued for many years, using a variety of methodologies.

One of the earliest methods was to isolate and analyze levels of the proteins in a cell, tissue, organ system, or even organisms both before and after the process of

interest. One method of analyzing multiple proteins in a sample is using 2-dimensional gel electrophoresis, wherein proteins can be, in principle, identified and quantified as individual bands, and ultimately reduced to a discrete signal. At present, 2-dimensional analysis only resolves approximately 15% of the proteins. In order to positively analyze those bands which are resolved, each band must be

approximately 15% of the proteins. In order to positively analyze those bands which are resolved, each band must be excised from the membrane and subjected to protein sequence analysis using Edman degradation. Unfortunately, most of the bands were present in quantities too small to obtain a reliable sequence, and many of those bands contained more

than one discrete protein. An additional difficulty is that many of the proteins were blocked at the amino-terminus, further complicating the sequencing process.

Analyzing differentiation at the gene transcription level has overcome many of these disadvantages and drawbacks, since the power of recombinant DNA technology allows amplification of signals containing very small 5 amounts of material. The most common method, called "hybridization subtraction," involves isolation of mRNA from the biological specimen before (B) and after (A) the developmental process of interest, transcribing one set of mRNA into cDNA, subtracting specimen B from specimen A 10 (mRNA from cDNA) by hybridization, and constructing a cDNA library from the non-hybridizing mRNA fraction. Many different groups have used this strategy successfully, and a variety of procedures have been published and improved upon using this same basic scheme. Nucl. Acids Res. 19, 15 7097-7104 (1991); Nucl. Acids Res. 18, 4833-42 (1990); · Nucl. Acids Res. 18, 2789-92 (1989); European J. Neuroscience 2, 1063-1073 (1990); Analytical Biochem. 187, 364-73 (1990); Genet. Annals Techn. Appl. 7, 64-70 (1990); . GATA 8(4), 129-33 (1991); Proc. Natl. Acad. Sci. USA 85, 1696-1700 (1988); Nucl. Acids Res. 19, 1954 (1991); Proc. Natl. Acad. Sci. USA 88, 1943-47 (1991); Nucl. Acids Res. 19, 6123-27 (1991); Proc. Natl. Acad. Sci. USA 85, 5738-42 (1988); Nucl. Acids Res. 16, 10937 (1988).

Although each of these techniques have particular 25 strengths and Weaknesses, there are still some limitations and undesirable aspects of these methods: First, the time and effort required to construct such libraries is quite large. Typically, a trained molecular biologist might expect construction and characterization of such a library 30 to require 3 to 6 months, depending on the level of skill, experience, and luck. Second, the resulting subtraction libraries are typically inferior to the libraries constructed by standard methodology. A typical conventional cDNA library should have a clone complexity of 35 at least 106 clones, and an average insert size of 1-3 kB. In contrast, subtracted libraries can have complexities of 10² or 10³ and average insert sizes of 0.2 kB. Therefore, there can be a significant loss of clone and sequence information associated with such libraries. Third, this

approach allows the researcher to capture only the genes induced in sp cimen A relative to specimen B, not vice-versa, nor does it easily allow comparison to a third specimen of interest (C). Fourth, this approach requires very large amounts (hundreds of micrograms) of "driver" mRNA (specimen B), which significantly limits the number and type of subtractions that are possible since many tissues and cells are very difficult to obtain in large quantities.

Fifth, the resolution of the subtraction is dependent 10 upon the physical properties of DNA: DNA or RNA: DNA hybridization. The ability of a given sequence to find a hybridization match is dependent on its unique CoT value. The CoT value is a function of the number of copies 15 (concentration) of the particular sequence, multiplied by the time of hybridization. It follows that for sequences which are abundant, hybridization events will occur very rapidly (low CoT value), while rare sequences will form duplexes at very high CoT values. CoT values which allow 20 such rare sequences to form duplexes and therefore be effectively selected are difficult to achieve in a convenient time frame. Therefore, hybridization subtraction is simply not a useful technique with which to study relative levels of rare mRNA species. Sixth, this 25 problem is further complicated by the fact that duplex formation is also dependent on the nucleotide base composition for a given sequence. Those sequences rich in G + C form stronger duplexes than those with high contents of A + T. Therefore, the former sequences will tend to be 30 removed selectively by hybridization subtraction. Seventh, it is possible that hybridization between nonexact matches can occur. When this happens, the expression of a homologous gene may "mask" expression of a gene of interest, artificially skewing the results for that 35 particular gene.

Matsubara and Okubo proposed using partial cDNA sequences to establish expression profiles of genes which could be used in functional analyses of the human genome.

Matsubara and Okubo warned against using random priming, as

it creates multiple unique DNA fragments from individual mRNAs and may thus skew the analysis of the number of particular mRNAs per library. They sequenced randomly selected members from a 3'-directed cDNA library and 5 established the frequency of appearance of the various ESTs. They proposed comparing lists of ESTs from various cell types to classify genes. Genes expressed in many different cell types were labeled housekeepers and those selectively expressed in certain cells were labeled cell-10 specific genes, even in the absence of the full sequence of the gene or the biological activity of the gene product.

The present invention avoids the drawbacks of the prior art by providing a method to quantify the relative abundance of multiple gene transcripts in a given

15 biological specimen by the use of high-throughput sequence-specific analysis of individual RNAs and/or their corresponding cDNAs.

The present invention offers several advantages over current protein discovery methods which attempt to isolate individual proteins based upon biological effects. The method of the instant invention provides for detailed diagnostic comparisons of cell profiles revealing numerous changes in the expression of individual transcripts.

The instant invention provides several advantages over

25 current subtraction methods including a more complex
library analysis (106 to 107 clones as compared to 103
clones) which allows identification of low abundance
messages as well as enabling the identification of messages
which either increase or decrease in abundance. These

30 large libraries are very routine to make in contrast to the
libraries of previous methods. In addition, homologues can
easily be distinguished with the method of the instant
invention.

This method is very convenient because it organizes a large quantity of data into a comprehensible, digestible format. The most significant differences are highlighted by electronic subtraction. In depth analyses are made more convenient.

The present invention provides s veral advantages over previous methods of electronic analysis of cDNA. The method is particularly powerful when more than 100 and preferably more than 1,000 gene transcripts are analyzed.

5 In such a case, new low-frequency transcripts are discovered and tissue typed.

High resolution analysis of gene expression can be used directly as a diagnostic profile or to identify disease-specific genes for the development of more classic diagnostic approaches.

This process is defined as gene transcript frequency analysis. The resulting quantitative analysis of the gene transcripts is defined as comparative gene transcript analysis.

3. SUMMARY OF THE INVENTION

15

The invention is a method of analyzing a specimen containing gene transcripts comprising the steps of (a) producing a library of biological sequences; (b) generating a set of transcript sequences, where each of the transcript sequences in said set is indicative of a different one of the biological sequences of the library; (c) processing the transcript sequences in a programmed computer (in which a database of reference transcript sequences indicative of reference sequences is stored), to generate an identified 25 sequence value for each of the transcript sequences, where each said identified sequence value is indicative of sequence annotation and a degree of match between one of the biological sequences of the library and at least one of the reference sequences; and (d) processing each said identified sequence value to generate final data values indicative of the number of times each identified sequence value is present in the library.

The invention also includes a method of comparing two specimens containing gene transcripts. The first specimen is processed as described above. The second specimen is used to produce a second library of biological sequences, which is used to generate a second set of transcript sequences, where each of the transcript sequences in the

second set is indicative of one of the biological sequences of the second library. Then the second set of transcript sequences is processed in a programmed computer to generate. a second set of identified sequence values, namely the 5 further identified sequence values, each of which is indicative of a sequence annotation and includes a degree of match between one of the biological sequences of the second library and at least one of the reference sequences. The further identified sequence values are processed to 10 generate further final data values indicative of the number of times each further identified sequence value is present in the second library. The final data values from the first specimen and the further identified sequence values from the second specimen are processed to generate ratios 15 of transcript sequences, which indicate the differences in the number of gene transcripts between the two specimens.

In a further embodiment, the method includes quantifying the relative abundance of mRNA in a biological specimen by (a) isolating a population of mRNA transcripts from a biological specimen; (b) identifying genes from which the mRNA was transcribed by a sequence-specific method; (c) determining the numbers of mRNA transcripts corresponding to each of the genes; and (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts.

Also disclosed is a method of producing a gene transcript image analysis by first obtaining a mixture of mRNA, from which cDNA copies are made. The cDNA is inserted into a suitable vector which is used to transfect suitable host strain cells which are plated out and permitted to grow into clones, each cone representing a unique mRNA. A representative population of clones transfected with cDNA is isolated. Each clone in the population is identified by a sequence-specific method which identifies the gene from which the unique mRNA was transcribed. The number of times each gene is identified to a clone is determined to evaluate gene transcript abundance. The genes and their abundances are listed in order of abundance to produce a gene transcript image.

In a further embodiment, the relative abundance of the gene transcripts in one cell type or tissue is compared with the relative abundance of gene transcript numbers in a second cell type or tissue in order to identify the differences and similarities.

In a further embodiment, the method includes a system for analyzing a library of biological sequences including a means for receiving a set of transcript sequences, where each of the transcript sequences is indicative of a 10 different one of the biological sequences of the library; and a means for processing the transcript sequences in a computer system in which a database of reference transcript sequences indicative of reference sequences is stored, wherein the computer is programmed with software for . 15 generating an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of a sequence annotation and the degree of match between a different one of the biological sequences of the library and at least one of the reference 20 sequences, and for processing each said identified sequence value to generate final data values indicative of the number of times each identified sequence value is present in the library.

In essence, the invention is a method and system for 25 quantifying the relative abundance of gene transcripts in a biological specimen. The invention provides a method for comparing the gene transcript image from two or more different biological specimens in order to distinguish between the two specimens and identify one or more genes which are differentially expressed between the two Thus, this gene transcript image and its specimens. comparison can be used as a diagnostic. One embodiment of the method generates high-throughput sequence-specific analysis of multiple RNAs or their corresponding cDNAs: a gene transcript image. Another embodiment of the method produces the gene transcript imaging analysis by the use of high-throughput cDNA sequence analysis. In addition, two or more gene transcript images can be compared and used to detect or diagnose a particular biological state, disease,

or condition which is correlated to the relative abundance of gene transcripts in a given cell or population of cells.

4. DESCRIPTION OF THE TABLES AND DRAWINGS 4.1. TABLES

5 <u>Table 1</u> presents a detailed explanation of the letter codes utilized in Tables 2-5.

Table 2 lists the one hundred most common gene transcripts. It is a partial list of isolates from the HUVEC cDNA library prepared and sequenced as described below. The left-hand column refers to the sequence's order of abundance in this table. The next column labeled "number" is the clone number of the first HUVEC sequence identification reference matching the sequence in the "entry" column number. Isolates that have not been sequenced are not present in Table 2. The next column, labeled "N", indicates the total number of cDNAs which have the same degree of match with the sequence of the reference transcript in the "entry" column.

The column labeled "entry" gives the NIH GENBANK locus
name, which corresponds to the library sequence numbers.
The "s" column indicates in a few cases the species of the
reference sequence. The code for column "s" is given in
Table 1. The column labeled "descriptor" provides a plain
English explanation of the identity of the sequence
corresponding to the NIH GENBANK locus name in the "entry"
column.

Table 3 is a comparison of the top fifteen most abundant gene transcripts in normal monocytes and activated macrophage cells.

Table 4 is a detailed summary of library subtraction analysis summary comparing the THP-1 and human macrophage cDNA sequences. In Table 4, the same code as in Table 2 is used. Additional columns are for "bgfreq" (abundance number in the subtractant library), "rfend" (abundance number in the target library) and "ratio" (the target abundance number divided by the subtractant abundance number). As is clear from perusal of the table, when the abundance number in the subtractant library is "0", the

target abundance number is divided by 0.05. This is a way of obtaining a result (not possible dividing by 0) and distinguishing the result from ratios of subtractant numbers of 1.

Table 5 is the computer program, written in source code, for generating gene transcript subtraction profiles.

Table 6 is a partial listing of database entries used in the electronic northern blot analysis as provided by the present invention.

10

35

5

4.2. BRIEF DESCRIPTION OF THE DRAWINGS

Figure 1 is a chart summarizing data collected and stored regarding the library construction portion of sequence preparation and analysis.

15 Figure 2 is a diagram representing the sequence of operations performed by "abundance sort" software in a class of preferred embodiments of the inventive method.

Figure 3 is a block diagram of a preferred embodiment of the system of the invention.

Figure 4 is a more detailed block diagram of the 20 bioinformatics process from new sequence (that has already been sequenced but not identified) to printout of the transcript imaging analysis and the provision of database subscriptions.

25 5. DETAILED DESCRIPTION OF THE INVENTION

The present invention provides a method to compare the relative abundance of gene transcripts in different biological specimens by the use of high-throughput sequence-specific analysis of individual RNAs or their 30 corresponding cDNAs (or alternatively, of data representing other biological sequences). This process is denoted herein as gene transcript imaging. The quantitative analysis of the relative abundance for a set of gene transcripts is denoted herein as "gene transcript image analysis" or "gene transcript frequency analysis". present invention allows one to obtain a profile for gene transcription in any given population of cells or tissue from any type of organism. The invention can be applied to

obtain a profile of a specimen consisting of a single cell (or clones of a single cell), or of many cells, or of tissue more complex than a single cell and containing multiple cell types, such as liver.

The invention has significant advantages in the fields of diagnostics, toxicology and pharmacology, to name a few. A highly sophisticated diagnostic test can be performed on the ill patient in whom a diagnosis has not been made. A biological specimen consisting of the patient's fluids or tissues is obtained, and the gene transcripts are isolated and expanded to the extent necessary to determine their identity. Optionally, the gene transcripts can be converted to cDNA. A sampling of the gene transcripts are subjected to sequence-specific analysis and quantified.

These gene transcript sequence abundances are compared against reference database sequence abundances including

against reference database sequence abundances are compared against reference database sequence abundances including normal data sets for diseased and healthy patients. The patient has the disease(s) with which the patient's data set most closely correlates.

20

For example, gene transcript frequency analysis can be used to differentiate normal cells or tissues from diseased cells or tissues, just as it highlights differences between normal monocytes and activated macrophages in Table 3.

In toxicology, a fundamental question is which tests

25 are most effective in predicting or detecting a toxic
effect. Gene transcript imaging provides highly detailed
information on the cell and tissue environment, some of
which would not be obvious in conventional, less detailed
screening methods. The gene transcript image is a more

30 powerful method to predict drug toxicity and efficacy.
Similar benefits accrue in the use of this tool in
pharmacology. The gene transcript image can be used
selectively to look at protein categories which are
expected to be affected, for example, enzymes which

35 detoxify toxins.

In an alternative embodiment, comparative gene transcript frequency analysis is used to differentiate between cancer cells which respond to anti-cancer agents and those which do not respond. Examples of anti-cancer

agents are tamoxifen, vincristine, vinblastine,
p dophyllotoxins, etoposide, tenisposide, cisplatin,
biologic response modifiers such as interfer n, I1-2, GMCSF, enzymes, h rmones and the like. This method also
provides a means for sorting the gene transcripts by
functional category. In the case of cancer cells,
transcription factors or other essential regulatory
molecules are very important categories to analyze across
different libraries.

In yet another embodiment, comparative gene transcript frequency analysis is used to differentiate between control liver cells and liver cells isolated from patients treated with experimental drugs like FIAU to distinguish between pathology caused by the underlying disease and that caused by the drug.

In yet another embodiment, comparative gene transcript frequency analysis is used to differentiate between brain tissue from patients treated and untreated with lithium.

In a further embodiment, comparative gene transcript 20 frequency analysis is used to differentiate between cyclosporin and FK506-treated cells and normal cells.

In a further embodiment, comparative gene transcript frequency analysis is used to differentiate between virally infected (including HIV-infected) human cells and uninfected human cells. Gene transcript frequency analysis is also used to rapidly survey gene transcripts in HIV-resistant, HIV-infected, and HIV-sensitive cells. Comparison of gene transcript abundance will indicate the success of treatment and/or new avenues to study.

In a further embodiment, comparative gene transcript frequency analysis is used to differentiate between bronchial lavage fluids from healthy and unhealthy patients with a variety of ailments.

30

In a further embodiment, comparative gene transcript
frequency analysis is used to differentiate between cell,
plant, microbial and animal mutants and wild-type species.
In addition, the transcript abundance program is adapted to
permit the scientist to evaluate the transcription of one
gene in many differ nt tissues. Such comparisons could

identify deletion mutants which do not produce a gene product and point mutants which produce a less abundant or therwise different message. Such mutations can affect basic biochemical and pharmacological processes, such as mineral nutrition and metabolism, and can be isolated by means known to those skilled in the art. Thus, crops with improved yields, pest resistance and other factors can be developed.

In a further embodiment, comparative gene transcript

frequency analysis is used for an interspecies comparative
analysis which would allow for the selection of better
pharmacologic animal models. In this embodiment, humans
and other animals (such as a mouse), or their cultured
cells are treated with a specific test agent. The relative

sequence abundance of each cDNA population is determined.

If the animal test system is a good model, homologous genes
in the animal cDNA population should change expression
similarly to those in human cells. If side effects are
detected with the drug, a detailed transcript abundance
analysis will be performed to survey gene transcript
changes. Models will then be evaluated by comparing basic
physiological changes.

In a further embodiment, comparative gene transcript frequency analysis is used in a clinical setting to give a 25 highly detailed gene transcript profile of a patient's cells or tissue (for example, a blood sample). In particular, gene transcript frequency analysis is used to give a high resolution gene expression profile of a diseased state or condition.

In the preferred embodiment, the method utilizes high-throughput cDNA sequencing to identify specific transcripts of interest. The generated cDNA and deduced amino acid sequences are then extensively compared with GENBANK and other sequence data banks as described below. The method offers several advantages over current protein discovery by two-dimensional gel methods which try to identify individual proteins involved in a particular biological effect. Here, detailed comparisons of profiles of activated and inactive cells reveal numerous changes in

30

the expression of individual transcripts. After it is determined if the sequence is an "exact" match, similar or a non-match, the sequence is entered into a database. Next, the numbers of copies of cDNA corresponding to each 5 gene are tabulated. Although this can be done slowly and arduously, if at all, by human hand from a printout of all entries, a computer program is a useful and rapid way to tabulate this information. The numbers of cDNA copies (optionally divided by the total number of sequences in the 10 data set) provides a picture of the relative abundance of transcripts for each corresponding gene. The list of represented genes can then be sorted by abundance in the cDNA population. A multitude of additional types of comparisons or dimensions are possible and are exemplified 15 below.

An alternate method of producing a gene transcript image includes the steps of obtaining a mixture of test mRNA and providing a representative array of unique probes whose sequences are complementary to at least some of the test mRNAs. Next, a fixed amount of the test mRNA is added to the arrayed probes. The test mRNA is incubated with the probes for a sufficient time to allow hybrids of the test mRNA and probes to form. The mRNA-probe hybrids are detected and the quantity determined. The hybrids are identified by their location in the probe array. The quantity of each hybrid is summed to give a population number. Each hybrid quantity is divided by the population number to provide a set of relative abundance data termed a gene transcript image analysis.

30

6. EXAMPLES

The examples below are provided to illustrate the subject invention. These examples are provided by way of illustration and are not included for the purpose of limiting the invention.

35 6.1. TISSUE SOURCES AND CELL LINES

For analysis with the computer program claimed herein, biological sequences can be obtained from virtually any

Most popular are tissues obtained from the human Tissues can be obtained from any organ of the body, any age donor, any abnormality or any immortalized cell Immortal cell lines may be pref rred in some 5 instances because of their purity of cell type; other tissue samples invariably include mixed cell types. special technique is available to take a single cell (for example, a brain cell) and harness the cellular machinery to grow up sufficient cDNA for sequencing by the techniques 10 and analysis described herein (cf. U.S. Patent Nos. 5,021,335 and 5,168,038, which are incorporated by The examples given herein utilized the reference). following immortalized cell lines: monocyte-like U-937 cells, activated macrophage-like THP-1 cells, induced 15 vascular endothelial cells (HUVEC cells) and mast cell-like HMC-1 cells.

The U-937 cell line is a human histiocytic lymphoma cell line with monocyte characteristics, established from malignant cells obtained from the pleural effusion of a 20 patient with diffuse histiocytic lymphoma (Sundstrom, C. and Nilsson, K. (1976) Int. J. Cancer 17:565). one of only a few human cell lines with the morphology, cytochemistry, surface receptors and monocyte-like characteristics of histiocytic cells. These cells can be 25 induced to terminal monocytic differentiation and will express new cell surface molecules when activated with supernatants from human mixed lymphocyte cultures. Upon this type of in vitro activation, the cells undergo morphological and functional changes, including 30 augmentation of antibody-dependent cellular cytotoxicity (ADCC) against erythroid and tumor target cells (one of the principal functions of macrophages). Activation of U-937 cells with phorbol 12-myristate 13-acetate (PMA) in vitro stimulates the production of several compounds, including 35 prostaglandins, leukotrienes and platelet-activating factor (PAF), which are potent inflammatory mediators. 937 is a cell line that is well suited for the identification and isolation of gene transcripts associated with normal monocytes.

The HUVEC cell line is a normal, homogeneous, well characterized, early passage endothelial cell culture from human umbilical vein (Cell Systems Corp., 12815 NE 124th Street, Kirkland, WA 98034). Only gene transcripts from induced, or treated, HUVEC cells were sequenced. One batch of 1 X 10⁸ cells was treated for 5 hours with 1 U/ml rIL-1b and 100 ng/ml E.coli lipopolysaccharide (LPS) endotoxin prior to harvesting. A separate batch of 2 X 10⁸ cells was treated at confluence with 4 U/ml TNF and 2 U/ml interferon-gamma (IFN-gamma) prior to harvesting.

THP-1 is a human leukemic cell line with distinct monocytic characteristics. This cell line was derived from the blood of a 1-year-old boy with acute monocytic leukemia (Tsuchiya, S. et al. (1980) Int. J. Cancer: 171-76). 15 following cytological and cytochemical criteria were used to determine the monocytic nature of the cell line: 1) the presence of alpha-naphthyl butyrate esterase activity which could be inhibited by sodium fluoride; 2) the production of lysozyme; 3) the phagocytosis of latex particles and 20 sensitized SRBC (sheep red blood cells); and 4) the ability of mitomycin C-treated THP-1 cells to activate Tlymphocytes following ConA (concanavalin A) treatment. Morphologically, the cytoplasm contained small azurophilic granules and the nucleus was indented and irregularly 25 shaped with deep folds. The cell line had Fc and C3b receptors, probably functioning in phagocytosis. cells treated with the tumor promoter 12-o-tetradecanoylphorbol-13 acetate (TPA) stop proliferating and differentiate into macrophage-like cells which mimic native 30 monocyte-derived macrophages in several respects. Morphologically, as the cells change shape, the nucleus becomes more irregular and additional phagocytic vacuoles appear in the cytoplasm. The differentiated THP-1 cells also exhibit an increased adherence to tissue culture 35 plastic.

HMC-1 cells (a human mast cell line) were established from the peripheral blood of a Mayo Clinic patient with mast cell leukemia (Leukemia Res. (1988) 12:345-55). The cultured cells looked similar to immature cloned murine

mast cells, contained histamine, and stained positively for chloroacetate esterase, amino caproate esterase, eosinophil major basic protein (MBP) and tryptase. The HMC-1 cells have, however, lost the ability to synthesize normal IgE receptors. HMC-1 cells also possess a 10;16 translocation, present in cells initially collected by leukophoresis from the patient and not an artifact of culturing. Thus, HMC-1 cells are a good model for mast cells.

6.2. CONSTRUCTION OF CDNA LIBRARIES

10 For inter-library comparisons, the libraries must be prepared in similar manners. Certain parameters appear to be particularly important to control. One such parameter is the method of isolating mRNA. It is important to use the same conditions to remove DNA and heterogeneous nuclear RNA from comparison libraries. Size fractionation of cDNA must be carefully controlled. The same vector preferably should be used for preparing libraries to be compared. At the very least, the same type of vector (e.g., unidirectional vector) should be used to assure a valid comparison. A unidirectional vector may be preferred in order to more easily analyze the output.

It is preferred to prime only with oligo dT
unidirectional primer in order to obtain one only clone per
mRNA transcript when obtaining cDNAs. However, it is

25 recognized that employing a mixture of oligo dT and random
primers can also be advantageous because such a mixture
results in more sequence diversity when gene discovery also
is a goal. Similar effects can be obtained with DR2
(Clontech) and HXLOX (US Biochemical) and also vectors from

30 Invitrogen and Novagen. These vectors have two
requirements. First, there must be primer sites for
commercially available primers such as T3 or M13 reverse
primers. Second, the vector must accept inserts up to 10
kB.

It also is important that the clones be randomly sampled, and that a significant population of clones is used. Data hav been generated with 5,000 clones; how ver, if very rare genes are to be obtained and/or their relativ

35

5

35

abundance determined, as many as 100,000 clones from a single library may need to be sampled. Size fractionation of cDNA also must be carefully controlled. Alternately, plagues can be selected, rather than clones.

Besides the Uni-ZAP™ vector system by Stratagene disclosed below, it is now believed that other similarly unidirectional vectors also can be used. For example, it is believed that such vectors include but are not limited to DR2 (Clontech), and HXLOX (U.S. Biochemical).

Preferably, the details of library construction (as shown in Figure 1) are collected and stored in a database for later retrieval relative to the sequences being compared. Fig. 1 shows important information regarding the library collaborator or cell or cDNA supplier,

15 pretreatment, biological source, culture, mRNA preparation and cDNA construction. Similarly detailed information about the other steps is beneficial in analyzing sequences and libraries in depth.

RNA must be harvested from cells and tissue samples
and cDNA libraries are subsequently constructed. cDNA
libraries can be constructed according to techniques known
in the art. (See, for example, Maniatis, T. et al. (1982)
Molecular Cloning, Cold Spring Harbor Laboratory, New
York). cDNA libraries may also be purchased. The U-937
cDNA library (catalog No. 937207) was obtained from
Stratagene, Inc., 11099 M. Torrey Pines Rd., La Jolla, CA
92037.

The THP-1 cDNA library was custom constructed by Stratagene from THP-1 cells cultured 48 hours with 100 nm TPA and 4 hours with 1 μ g/ml LPS. The human mast cell HMC-1 cDNA library was also custom constructed by Stratagene from cultured HMC-1 cells. The HUVEC cDNA library was custom constructed by Stratagene from two batches of induced HUVEC cells which were separately processed.

Essentially, all the libraries were prepared in the same manner. First, poly(A+)RNA (mRNA) was purified. For the U-937 and HMC-1 RNA, cDNA synthesis was only primed with oligo dT. For the THP-1 and HUVEC RNA, cDNA synthesis was primed separately with both oligo dT and random

hexamers, and the two cDNA libraries were treated separately. Synthetic adaptor oligonucleotides were ligated onto cDNA ends enabling its insertion into th Uni-. Zap™ vector system (Stratagene), allowing high efficiency 5 unidirectional (sense orientation) lambda library construction and the convenience of a plasmid system with blue-white color selection to detect clones with cDNA insertions. Finally, the two libraries were combined into a single library by mixing equal numbers of bacteriophage.

The libraries can be screened with either DNA probes or antibody probes and the pBluescript® phagemid (Stratagene) can be rapidly excised in vivo. The phagemid allows the use of a plasmid system for easy insert characterization, sequencing, site-directed mutagenesis, 15 the creation of unidirectional deletions and expression of fusion proteins. The custom-constructed library phage particles were infected into E. coli host strain XL1-Blue® (Stratagene), which has a high transformation efficiency, increasing the probability of obtaining rare, under-20 represented clones in the cDNA library.

10

6.3. ISOLATION OF CDNA CLONES

The phagemid forms of individual cDNA clones were obtained by the in vivo excision process, in which the host bacterial strain was coinfected with both the lambda library phage and an f1 helper phage. Proteins derived from both the library-containing phage and the helper phage nicked the lambda DNA, initiated new DNA synthesis from defined sequences on the lambda target DNA and created a smaller, single stranded circular phagemid DNA molecule 30 that included all DNA sequences of the pBluescript® plasmid and the cDNA insert. The phagemid DNA was secreted from the cells and purified, then used to re-infect fresh host cells, where the double stranded phagemid DNA was produced. Because the phagemid carries the gene for beta-lactamase, 35 the newly-transformed bacteria are selected on medium containing ampicillin.

Phagemid DNA was purified using the Magic Minipreps™ DNA Purification System (Promega catalogue #A7100. Promega

Corp., 2800 Woods Hollow Rd., Madison, WI 53711). This small-scale process provides a simple and reliable method for lysing the bacterial cells and rapidly isolating purified phagemid DNA using a proprietary DNA-binding resin. The DNA was eluted from the purification resin already prepared for DNA sequencing and other analytical manipulations.

Phagemid DNA was also purified using the QIAwell-8
Plasmid Purification System from QIAGEN® DNA Purification

System (QIAGEN Inc., 9259 Eton Ave., Chattsworth, CA
91311). This product line provides a convenient, rapid and
reliable high-throughput method for lysing the bacterial
cells and isolating highly purified phagemid DNA using
QIAGEN anion-exchange resin particles with EMPORE™ membrane

technology from 3M in a multiwell format. The DNA was
eluted from the purification resin already prepared for DNA
sequencing and other analytical manipulations.

An alternate method of purifying phagemid has recently become available. It utilizes the Miniprep Kit (Catalog 20 No. 77468, available from Advanced Genetic Technologies Corp., 19212 Orbit Drive, Gaithersburg, Maryland). kit is in the 96-well format and provides enough reagents for 960 purifications. Each kit is provided with a recommended protocol, which has been employed except for 25 the following changes. First, the 96 wells are each filled with only 1 ml of sterile terrific broth with carbenicillin at 25 mg/L and glycerol at 0.4%. After the wells are inoculated, the bacteria are cultured for 24 hours and lysed with 60 μ l of lysis buffer. A centrifugation step 30 (2900 rpm for 5 minutes) is performed before the contents of the block are added to the primary filter plate. optional step of adding isopropanol to TRIS buffer is not routinely performed. After the last step in the protocol, samples are transferred to a Beckman 96-well block for 35 storage.

Another new DNA purification system is the WIZARD™ product line which is available from Promega (catalog No. A7071) and may be adaptable to the 96-w ll format.

6.4. SEQUENCING OF CONA CLONES

The cDNA inserts from random isolates of the U-937 and THP-1 libraries were sequenced in part. Methods for DNA sequencing are well known in the art. Conventional 5 enzymatic methods employ DNA polymerase Klenow fragment, Sequenase™ or Taq polymerase to extend DNA chains from an oligonucleotide primer annealed to the DNA template of interest. Methods have been developed for the use of both single- and double-stranded templates. The chain 10 termination reaction products are usually electrophoresed on urea-acrylamide gels and are detected either by autoradiography (for radionuclide-labeled precursors) or by fluorescence (for fluorescent-labeled precursors). Recent improvements in mechanized reaction preparation, sequencing 15 and analysis using the fluorescent detection method have permitted expansion in the number of sequences that can be determined per day (such as the Applied Biosystems 373 and 377 DNA sequencer, Catalyst 800). Currently with the system as described, read lengths range from 250 to 400 20 bases and are clone dependent. Read length also varies with the length of time the gel is run. In general, the shorter runs tend to truncate the sequence. A minimum of only about 25 to 50 bases is necessary to establish the identification and degree of homology of the sequence. 25 Gene transcript imaging can be used with any sequencespecific method, including, but not limited to hybridization, mass spectroscopy, capillary electrophoresis and 505 gel electrophoresis.

6.5. HOMOLOGY SEARCHING OF CDNA CLONE AND DEDUCED PROTEIN (and Subsequent Steps)

30

Using the nucleotide sequences derived from the cDNA clones as query sequences (sequences of a Sequence Listing), databases containing previously identified sequences are searched for areas of homology (similarity).

35 Examples of such databases include Genbank and EMBL. We next describe examples of two homology search algorithms that can be used, and then describe the subsequent computer-implemented steps to be performed in accordance with preferred embodiments of the invention.

In the following description of the computerimplemented steps of the invention, the word "library"
denotes a set (or population) of biological specimen
nucleic acid sequences. A "library" can consist of cDNA
sequences, RNA sequences, or the like, which characterize a
biological specimen. The biological specimen can consist
of cells of a single human cell type (or can be any of the
other above-mentioned types of specimens). We contemplate
that the sequences in a library have been determined so as
to accurately represent or characterize a biological
specimen (for example, they can consist of representative
cDNA sequences from clones of RNA taken from a single human
cell).

In the following description of the computerimplemented steps of the invention, the expression
"database" denotes a set of stored data which represent a
collection of sequences, which in turn represent a
collection of biological reference materials. For example,
a database can consist of data representing many stored
cDNA sequences which are in turn representative of human
cells infected with various viruses, cells of humans of
various ages, cells from different mammalian species, and
so on.

In preferred embodiments, the invention employs a computer programmed with software (to be described) for performing the following steps:

- (a) processing data indicative of a library of cDNA sequences (generated as a result of high-throughput cDNA sequencing or other method) to determine whether each sequence in the library matches a DNA sequence of a reference database of DNA sequences (and if so, identifying the reference database entry which matches the sequence and indicating the degree of match between the reference sequence and the library sequence) and assigning an identified sequence value based on the sequence annotation and degree of match to each of the sequences in the library;
 - (b) for some or all entries of the database, tabulating the number of matching identified sequence

values in the library (Although this can be done by human hand from a printout of all entries, we prefer to perform this step using computer software to be described below.), thereby generating a set of final data values or "abundance 5 numbers"; and

(c) if the libraries are different sizes, dividing each abundance number by the total number of sequences in the library, to obtain a relative abundance number for each identified sequence value (i.e., a relative abundance of 10 each gene transcript).

The list of identified sequence values (or genes corresponding thereto) can then be sorted by abundance in the cDNA population. A multitude of additional types of comparisons or dimensions are possible.

15

For example (to be described below in greater detail), steps (a) and (b) can be repeated for two different libraries (sometimes referred to as a "target" library and a "subtractant" library). Then, for each identified sequence value (or gene transcript), a "ratio" value is 20 obtained by dividing the abundance number (for that identified sequence value) for the target library, by the abundance number (for that identified sequence value) for the subtractant library.

In fact, subtraction may be carried out on multiple 25 libraries. It is possible to add the transcripts from several libraries (for example, three) and then to divide them by another set of transcripts from multiple libraries (again, for example, three). Notation for this operation may be abbreviated as (A+B+C) / (D+E+F), where the capital 30 letters each indicate an entire library. Optionally the abundance numbers of transcripts in the summed libraries may be divided by the total sample size before subtraction.

Unlike standard hybridization technology which permits a single subtraction of two libraries, once one has 35 processed a set or library transcript sequences and stored them in the computer, any number of subtractions can be performed on the library. For example, by this method, ratio values can be obtained by dividing relative abundance

values in a first library by corresponding values in a second library and vice versa.

In variations on step (a), the library consists of nucleotide sequences derived from cDNA clones. Examples of 5 databases which can be searched for areas of homology (similarity) in step (a) include the commercially available databases known as Genbank (NIH) EMBL (European Molecular Biology Labs, Germany), and GENESEQ (Intelligenetics, Mountain View, California).

One homology search algorithm which can be used to 10. implement step (a) is the algorithm described in the paper by D.J. Lipman and W.R. Pearson, entitled "Rapid and Sensitive Protein Similarity Searches, "Science, 227:1435 (1985). In this algorithm, the homologous regions are 15 searched in a two-step manner. In the first step, the highest homologous regions are determined by calculating a matching score using a homology score table. The parameter "Ktup" is used in this step to establish the minimum window size to be shifted for comparing two sequences. Ktup also 20 sets the number of bases that must match to extract the highest homologous region among the sequences. step, no insertions or deletions are applied and the homology is displayed as an initial (INIT) value.

In the second step, the homologous regions are aligned 25 to obtain the highest matching score by inserting a gap in order to add a probable deleted portion. The matching score obtained in the first step is recalculated using the homology score Table and the insertion score Table to an optimized (OPT) value in the final output.

DNA homologies between two sequences can be examined graphically using the Harr method of constructing dot matrix homology plots (Needleman, S.B. and Wunsch, C.O., J. Mom. Biol 48:443 (1970)). This method produces a two-dimensional plot which can be useful in determining 35 regions of homology versus regions of repetition.

30

However, in a class of preferred embodiments, step (a) is implemented by processing the library data in the commercially available computer program known as the INHERIT 670 Sequence Analysis System, available from

Applied Biosystems Inc. (Foster City, California), including the software known as the Factura software (also available from Applied Biosystems Inc.). The Factura program preprocesses each library sequence to "edit out" 5 portions thereof which are not likely to be of interest, such as the vector used to prepare the library. Additional sequences which can be edited out or masked (ignored by the search tools) include but are not limited to the polyA tail and repetitive GAG and CCC sequences. A low-end search 10 program can be written to mask out such "low-information" sequences, or programs such as BLAST can ignore the lowinformation sequences.

In the algorithm implemented by the INHERIT 670 Sequence Analysis System, the Pattern Specification 15 Language (developed by TRW Inc.) is used to determine regions of homology. "There are three parameters that determine how INHERIT analysis runs sequence comparisons: window size, window offset and error tolerance. Window size specifies the length of the segments into which the 20 query sequence is subdivided. Window offset specifies where to start the next segment [to be compared], counting from the beginning of the previous segment. tolerance specifies the total number of insertions, deletions and/or substitutions that are tolerated over the specified word length. Error tolerance may be set to any integer between 0 and 6. The default settings are window tolerance=20, window offset=10 and error tolerance=3." INHERIT Analysis Users Manual, pp.2-15. Version 1.0, Applied Biosystems, Inc., October 1991.

Using a combination of these three parameters, a database (such as a DNA database) can be searched for sequences containing regions of homology and the appropriate sequences are scored with an initial value. Subsequently, these homologous regions are examined using 35 dot matrix homology plots to determine regions of homology versus regions of repetition. Smith-Waterman alignments can be used to display the results of the homology search. The INHERIT software can be executed by a Sun computer system programmed with the UNIX operating system.

30

Search alternatives to INHERIT include the BLAST program, GCG (available from the Gen tics Computer Group, WI) and the Dasher program (Temple Smith, Boston University, Boston, MA). Nucleotide sequences can be 5 searched against Genbank, EMBL or custom databases such as GENESEQ (available from Intelligenetics, Mountain View, CA) or other databases for genes. In addition, we have searched some sequences against our own in-house database.

In preferred embodiments, the transcript sequences are 10 analyzed by the INHERIT software for best conformance with a reference gene transcript to assign a sequence identifier and assigned the degree of homology, which together are the identified sequence value and are input into, and further processed by, a Macintosh personal computer (available from 15 Apple) programmed with an "abundance sort and subtraction analysis" computer program (to be described below).

Prior to the abundance sort and subtraction analysis program (also denoted as the "abundance sort" program), identified sequences from the cDNA clones are assigned 20 value (according to the parameters given above) by degree of match according to the following categories: "exact" matches (regions with a high degree of identity), homologous human matches (regions of high similarity, but not "exact" matches), homologous non-human matches (regions 25 of high similarity present in species other than human), or non matches (no significant regions of homology to previously identified nucleotide sequences stored in the form of the database). Alternately, the degree of match can be a numeric value as described below.

With reference again to the step of identifying matches between reference sequences and database entries, protein and peptide sequences can be deduced from the nucleic acid sequences. Using the deduced polypeptide sequence, the match identification can be performed in a 35 manner analogous to that done with cDNA sequences. protein sequence is used as a query sequence and compared to the previously identified sequences contained in a database such as the Swiss/Prot, PIR and the NBRF Protein database to find homologous proteins. These proteins are

30

initially scored for homology using a homology score Table (Orcutt, B.C. and Dayoff, M.O. Scoring Matrices, PIR Report MAT - 0285 (February 1985)) resulting in an INIT score. The homologous regions are aligned to obtain the highest matching scores by inserting a gap which adds a probable deleted portion. The matching score is recalculated using the homology score Table and the insertion score Table resulting in an optimized (OPT) score. Even in the absence of knowledge of the proper reading frame of an isolated sequence, the above-described protein homology search may be performed by searching all 3 reading frames.

Peptide and protein sequence homologies can also be ascertained using the INHERIT 670 Sequence Analysis System 15 in an analogous way to that used in DNA sequence homologies. Pattern Specification Language and parameter windows are used to search protein databases for sequences containing regions of homology which are scored with an initial value. Subsequent display in a dot-matrix homology 20 plot shows regions of homology versus regions of repetition. Additional search tools that are available to use on pattern search databases include PLsearch Blocks (available from Henikoff & Henikoff, University of Washington, Seattle), Dasher and GCG. Pattern search databases include, but are not limited to, Protein Blocks (available from Henikoff & Henikoff, University of Washington, Seattle), Brookhaven Protein (available from the Brookhaven National Laboratory, Brookhaven, MA), PROSITE (available from Amos Bairoch, University of Geneva, Switzerland), ProDom (available from Temple Smith, Boston University), and PROTEIN MOTIF FINGERPRINT (available from University of Leeds, United Kingdom).

The ABI Assembler application software, part of the INHERIT DNA analysis system (available from Applied Biosystems, Inc., Foster City, CA), can be employed to create and manage sequence assembly projects by assembling data from selected sequence fragments into a larger sequence. The Assembler software combines two advanced computer technologies which maximize the ability to

assemble sequenced DNA fragments into Assemblages, a special grouping of data where the relationships between sequences are shown by graphic overlap, alignment and statistical views. The process is based on the

Meyers-Kececioglu model of fragment assembly (INHERIT**
Assembler User's Manual, Applied Biosystems, Inc., Foster City, CA), and uses graph theory as the foundation of a

very rigorous multiple sequence alignment engine for assembling DNA sequence fragments. Other assembly programs that can be used include MEGALIGN (available from DNASTAR Inc., Madison, WI), Dasher and STADEN (available from Roger Staden, Cambridge, England).

Next, with reference to Fig. 2, we describe in more detail the "abundance sort" program which implements above15 mentioned "step (b)" to tabulate the number of sequences of the library which match each database entry (the "abundance number" for each database entry).

Fig. 2 is a flow chart of a preferred embodiment of the abundance sort program. A source code listing of this embodiment of the abundance sort program is set forth in Table 5. In the Table 5 implementation, the abundance sort program is written using the FoxBASE programming language commercially available from Microsoft Corporation. Although FoxBASE was the program chosen for the first iteration of this technology, it should not be considered limiting. Many other programming languages, Sybase being a particularly desirable alternative, can also be used, as will be obvious to one with ordinary skill in the art. The subroutine names specified in Fig. 2 correspond to

With reference again to Fig. 2, the "Identified Sequences" are transcript sequences representing each sequence of the library and a corresponding identification of the database entry (if any) which it matches. In other words, the "Identified Sequences" are transcript sequences representing the output of above-discussed "step (a)."

30 subroutines listed in Table 5.

Fig. 3 is a block diagram of a system for implementing the invention. The Fig. 3 system includes library generation unit 2 which generates a library and asserts an

output stream of transcript sequences indicative of the biological sequences comprising the library. Programmed processor 4 receives the data stream output from unit 2 and processes this data in accordance with above-discussed step (a) to generate the Identified Sequences. Processor 4 can be a processor programmed with the commercially available computer program known as the INHERIT 670 Sequence Analysis System and the commercially available computer program known as the Factura program (both available from Applied Biosystems Inc.) and with the UNIX operating system.

Still with reference to Fig. 3, the Identified Sequences are loaded into processor 6 which is programmed with the abundance sort program. Processor 6 generates the Final Transcript sequences indicated in both Figs. 2 and 3. Fig. 4 shows a more detailed block diagram of a planned relational computer system, including various searching techniques which can be implemented, along with an assortment of databases to guery against.

With reference to Fig. 2, the abundance sort program 20 first performs an operation known as "Tempnum" on the Identified Sequences, to discard all of the Identified Sequences except those which match database entries of selected types. For example, the Tempnum process can 25 select Identified Sequences which represent matches of the following types with database entries (see above for definition): "exact" matches, human "homologous" matches, "other species" matches representing genes present in species other than human), "no" matches (no significant 30 regions of homology with database entries representing previously identified nucleotide sequences), "I" matches (Incyte for not previously known DNA sequences), or "X" matches (matches ESTs in reference database). eliminates the U, S, M, V, A, R and D sequence (see Table 1 35 for definitions).

The identified sequence values selected during the "Tempnum" process then undergo a further selection (weeding out) operation known as "Tempr d." This operation can, for

example, discard all identified sequence values representing matches with selected database entries.

The identified sequence values selected during the "Tempred" process are then classified according to library, during the "Tempdesig" operation. It is contemplated that the "Identified Sequences" can represent sequences from a single library, or from two or more libraries.

Consider first the case that the identified sequence values represent sequences from a single library. 10 case, all the identified sequence values determined during "Tempred" undergo sorting in the "Templib" operation, further sorting in the "Libsort" operation, and finally additional sorting in the "Temptarsort" operation. For example, these three sorting operations can sort the 15 identified sequences in order of decreasing "abundance number" (to generate a list of decreasing abundance numbers, each abundance number corresponding to a unique identified sequence entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list 20 corresponding to database entries of a selected type) with redundancies eliminated from each sorted list. case, the operation identified as "Cruncher" can be bypassed, so that the "Final Data" values are the organized transcript sequences produced during the "Temptarsort" 25 operation.

We next consider the case that the transcript sequences produced during the "Tempred" operation represent sequences from two libraries (which we will denote the "target" library and the "subtractant" library). For 30 example, the target library may consist of cDNA sequences from clones of a diseased cell, while the subtractant library may consist of cDNA sequences from clones of the diseased cell after treatment by exposure to a drug. For another example, the target library may consist of cDNA sequences from clones of a cell type from a young human, while the subtractant library may consist of cDNA sequences from clones of the same cell type from the same human at different ages.

In this case, the "Tempdesig" operation routes all transcript sequences representing the target library for processing in accordance with "Templib" (and then "Libsort" and "Temptarsort"), and routes all transcript sequences 5 representing the subtractant library for processing in accordance with "Tempsub" (and then "Subsort" and "Tempsubsort"). For example, the consecutive "Templib," "Libsort," and "Temptarsort" sorting operations sort identified sequences from the target library in order of 10 decreasing abundance number (to generate a list of decreasing abundance numbers, each abundance number corresponding to a database entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list corresponding to database entries of a selected 15 type) with redundancies eliminated from each sorted list. The consecutive "Tempsub," "Subsort," and "Tempsubsort" sorting operations sort identified sequences from the subtractant library in order of decreasing abundance number (to generate a list of decreasing abundance numbers, each 20 abundance number corresponding to a database entry, or several lists of decreasing abundance numbers, with the abundance numbers in each list corresponding to database entries of a selected type) with redundancies eliminated from each sorted list.

operation typically represent sorted lists from which a histogram could be generated in which position along one (e.g., horizontal) axis indicates abundance number (of target library sequences), and position along another

(e.g., vertical) axis indicates identified sequence value (e.g., human or non-human gene type). Similarly, the transcript sequences output from the "Tempsubsort" operation typically represent sorted lists from which a histogram could be generated in which position along one

(e.g., horizontal) axis indicates abundance number (of subtractant library sequences), and position along another (e.g., vertical) axis indicates identified sequence value (e.g., human or non-human g ne type).

The transcript sequences (sorted lists) output from the Tempsubsort and Temptarsort sorting operations are combined during the operation identified as "Cruncher." The "Cruncher" process identifies pairs of corresponding 5 target and subtractant abundance numbers (both representing the same identified sequence value), and divides one by the other to generate a "ratio" value for each pair of corresponding abundance numbers, and then sorts the ratio values in order of decreasing ratio value. The data output 10 from the "Cruncher" operation (the Final Transcript sequence in Fig. 2) is typically a sorted list from Which a histogram could be generated in which position along one axis indicates the size of a ratio of abundance numbers (for corresponding identified sequence values from target 15 and subtractant libraries) and position along another axis indicates identified sequence value (e.g., gene type).

Preferably, prior to obtaining a ratio between the two library abundance values, the Cruncher operation also divides each ratio value by the total number of sequences in one or both of the target and subtractant libraries. The resulting lists of "relative" ratio values generated by the Cruncher operation are useful for many medical, scientific, and industrial applications. Also preferably, the output of the Cruncher operation is a set of lists, each list representing a sequence of decreasing ratio values for a different selected subset (e.g. protein family) of database entries.

In one example, the abundance sort program of the invention tabulates for a library the numbers of mRNA

30 transcripts corresponding to each gene identified in a database. These numbers are divided by the total number of clones sampled. The results of the division reflect the relative abundance of the mRNA transcripts in the cell type or tissue from which they were obtained. Obtaining this final data set is referred to herein as "gene transcript image analysis." The resulting subtracted data show exactly what proteins and genes are upregulated and downregulated in highly detailed complexity.

6.6. HUVEC CDNA LIBRARY

Table 2 is an abundance table listing the various gene transcripts in an induced HUVEC library. The transcripts are listed in order of decreasing abundance. This

5 computerized sorting simplifies analysis of the tissue and speeds identification of significant new proteins which are specific to this cell type. This type of endothelial cell lines tissues of the cardiovascular system, and the more that is known about its composition, particularly in

10 response to activation, the more choices of protein targets become available to affect in treating disorders of this tissue, such as the highly prevalent atherosclerosis.

6.7. MONOCYTE-CELL AND MAST-CELL CDNA LIBRARIES

Tables 3 and 4 show truncated comparisons of two 15 libraries. In Tables 3 and 4 the "normal monocytes" are the HMC-1 cells, and the "activated macrophages" are the THP-1 cells pretreated with PMA and activated with LPS. Table 3 lists in descending order of abundance the most abundant gene transcripts for both cell types. With only 15 gene transcripts from each cell type, this table permits 20 quick, qualitative comparison of the most common transcripts. This abundance sort, with its convenient side-by-side display, provides an immediately useful research tool. In this example, this research tool discloses that 1) only one of the top 15 activated macrophage transcripts is found in the top 15 normal monocyte gene transcripts (poly A binding protein); and 2) a new gene transcript (previously unreported in other databases) is relatively highly represented in activated 30 macrophages but is not similarly prominent in normal macrophages. Such a research tool provides researchers with a short-cut to new proteins, such as receptors, cellsurface and intracellular signalling molecules, which can serve as drug targets in commercial drug screening 35 programs. Such a tool could save considerable time over that consumed by a hit and miss discovery program aimed at identifying important proteins in and around cells, because those proteins carrying out everyday cellular functions and

represented as steady state mRNA are quickly eliminated from further characterization.

This illustrates how the g ne transcript profiles change with altered cellular function. Those skilled in 5 the art know that the biochemical composition of cells also changes with other functional changes such as cancer, including cancer's various stages, and exposure to toxicity. A gene transcript subtraction profile such as in Table 3 is useful as a first screening tool for such gene expression and protein studies.

6.8. SUBTRACTION ANALYSIS OF NORMAL MONOCYTE-CELL AND ACTIVATED MONOCYTE CELL CDNA LIBRARIES

Once the cDNA data are in the computer, the computer program as disclosed in Table 5 was used to obtain ratios 15 of all the gene transcripts in the two libraries discussed in Example 6.7, and the gene transcripts were sorted by the descending values of their ratios. If a gene transcript is not represented in one library, that gene transcript's abundance is unknown but appears to be less than 1. As an 20 approximation -- and to obtain a ratio, which would not be possible if the unrepresented gene were given an abundance of zero -- genes which are represented in only one of the two libraries are assigned an abundance of 1/2. Using 1/2 for unrepresented clones increases the relative importance 25 of "turned-on" and "turned-off" genes, whose products would be drug candidates. The resulting print-out is called a subtraction table and is an extremely valuable screening method, as is shown by the following data.

Table 4 is a subtraction table, in which the normal
monocyte library was electronically "subtracted" from the
activated macrophage library. This table highlights most
effectively the changes in abundance of the gene
transcripts by activation of macrophages. Even among the
first 20 gene transcripts listed, there are several unknown
gene transcripts. Thus, electronic subtraction is a useful
tool with which to assist researchers in identifying much
more quickly the basic biochemical changes between two cell
types. Such a tool can save universities and
pharmaceutical companies which spend billions of dollars on

research valuable time and laboratory resources at the early discovery stage and can speed up the drug development cycle, which in turn permits researchers to set up drug screening programs much earlier. Thus, this research tool provides a way to get new drugs to the public faster and more economically.

Also, such a subtraction table can be obtained for patient diagnosis. An individual patient sample (such as monocytes obtained from a biopsy or blood sample) can be compared with data provided herein to diagnose conditions associated with macrophage activation.

Table 4 uncovered many new gene transcripts (labeled Incyte clones). Note that many genes are turned on in the activated macrophage (i.e., the monocyte had a 0 in the bgfreq column). This screening method is superior to other screening techniques, such as the western blot, which are incapable of uncovering such a multitude of discrete new gene transcripts.

The subtraction-screening technique has also uncovered 20 a high number of cancer gene transcripts (oncogenes rho, ETS2, rab-2 ras, YPT1-related, and acute myeloid leukemia mRNA) in the activated macrophage. These transcripts may be attributed to the use of immortalized cell lines and are inherently interesting for that reason. This screening 25 technique offers a detailed picture of upregulated transcripts including oncogenes, which helps explain why anti-cancer drugs interfere with the patient's immunity mediated by activated macrophages. Armed with knowledge gained from this screening method, those skilled in the art 30 can set up more targeted, more effective drug screening programs to identify drugs which are differentially effective against 1) both relevant cancers and activated macrophage conditions with the same gene transcript profile; 2) cancer alone; and 3) activated macrophage 35 conditions.

Smooth muscle senescent protein (22 kd) was upregulated in the activated macrophage, which indicates that it is a candidate to block in controlling inflammation.

PCT/US95/01160 W 95/20681

SUBTRACTION ANALYSIS OF NORMAL LIVER CELLS AND HEPATITIS INFECTED LIVER CELL CDNA LIBRARIES

In this example, rats are exposed to hepatitis virus ' and maintained in the colony until they show definite signs 5 of hepatitis. Of the rats diagnosed with hepatitis, one half of the rats are treated with a new anti-hepatitis agent (AHA). Liver samples are obtained from all rats before exposure to the hepatitis virus and at the end of AHA treatment or no treatment. In addition, liver samples 10 can be obtained from rats with hepatitis just prior to AHA treatment.

The liver tissue is treated as described in Examples 6.2 and 6.3 to obtain mRNA and subsequently to sequence cDNA. The cDNA from each sample are processed and analyzed 15 for abundance according to the computer program in Table 5. The resulting gene transcript images of the cDNA provide detailed pictures of the baseline (control) for each animal and of the infected and/or treated state of the animals. cDNA data for a group of samples can be combined into a 20 group summary gene transcript profile for all control samples, all samples from infected rats and all samples from AHA-treated rats.

Subtractions are performed between appropriate individual libraries and the grouped libraries. individual animals, control and post-study samples can be subtracted. Also, if samples are obtained before and after AHA treatment, that data from individual animals and treatment groups can be subtracted. In addition, the data for all control samples can be pooled and averaged. control average can be subtracted from averages of both post-study AHA and post-study non-AHA cDNA samples. pre- and post-treatment samples are available, pre- and post-treatment samples can be compared individually (or electronically averaged) and subtracted.

These subtraction tables are used in two general ways. First, the differences are analyzed for gene transcripts which are associated with continuing hepatic deterioration or healing. The subtraction tables are tools to isolate the effects of the drug treatment from the underlying basic 40 pathology of hepatitis. Because hepatitis affects many

35

parameters, additional liver toxicity has been difficult to detect with only blood tests for the usual enzymes. The gene transcript profile and subtraction provides a much more complex biochemical picture which researchers have needed to analyze such difficult problems.

Second, the subtraction tables provide a tool for identifying clinical markers, individual proteins or other biochemical determinants which are used to predict and/or evaluate a clinical endpoint, such as disease, improvement 10 due to the drug, and even additional pathology due to the drug. The subtraction tables specifically highlight genes which are turned on or off. Thus, the subtraction tables provide a first screen for a set of gene transcript candidates for use as clinical markers. Subsequently, 15 electronic subtractions of additional cell and tissue libraries reveal which of the potential markers are in fact found in different cell and tissue libraries. Candidate gene transcripts found in additional libraries are removed from the set of potential clinical markers. Then, tests of 20 blood or other relevant samples which are known to lack and have the relevant condition are compared to validate the selection of the clinical marker. In this method, the particular physiologic function of the protein transcript need not be determined to qualify the gene transcript as a clinical marker.

6.10. ELECTRONIC NORTHERN BLOT

One limitation of electronic subtraction is that it is difficult to compare more than a pair of images at once.

Once particular individual gene products are identified as relevant to further study (via electronic subtraction or other methods), it is useful to study the expression of single genes in a multitude of different tissues. In the lab, the technique of "Northern" blot hybridization is used for this purpose. In this technique, a single cDNA, or a probe corresponding thereto, is labeled and then hybridized against a blot containing RNA samples prepared from a multitude of tissues or cell types. Upon autoradiography,

the pattern of expression of that particular gene, one at a time, can be quantitated in all the included samples.

In contrast, a further embodiment of this invention is the computerized form of this process, termed here

5 "electronic northern blot." In this variation, a single gene is queried for expression against a multitude of prepared and sequenced libraries present within the database. In this way, the pattern of expression of any single candidate gene can be examined instantaneously and effortlessly. More candidate genes can thus be scanned, leading to more frequent and fruitfully relevant discoveries. The computer program included as Table 5 includes a program for performing this function, and Table 6 is a partial listing of entries of the database used in the electronic northern blot analysis.

6.11. PHASE I CLINICAL TRIALS

Based on the establishment of safety and effectiveness in the above animal tests, Phase I clinical tests are undertaken. Normal patients are subjected to the usual 20 preliminary clinical laboratory tests. In addition, appropriate specimens are taken and subjected to gene transcript analysis. Additional patient specimens are taken at predetermined intervals during the test. The specimens are subjected to gene transcript analysis as described above. In addition, the gene transcript changes noted in the earlier rat toxicity study are carefully evaluated as clinical markers in the followed patients. Changes in the gene transcript analyses are evaluated as indicators of toxicity by correlation with clinical signs and symptoms and other laboratory results. In addition, subtraction is performed on individual patient specimens and on averaged patient specimens. The subtraction analysis highlights any toxicological changes in the treated patients. This is a highly refined determinant of 35 toxicity. The subtraction method also annotates clinical markers. Further subgroups can be analyzed by subtraction analysis, including, for example, 1) segregation by

occurrence and type of adverse effect; and 2) segregation by dosage.

6.12. GENE TRANSCRIPT IMAGING ANALYSIS IN CLINICAL STUDIES

A gene transcript imaging analysis (or multiple gene
transcript imaging analyses) is a useful tool in other
clinical studies. For example, the differences in gene
transcript imaging analyses before and after treatment can
be assessed for patients on placebo and drug treatment.
This method also effectively screens for clinical markers
to follow in clinical use of the drug.

6.13. COMPARATIVE GENE TRANSCRIPT ANALYSIS BETWEEN SPECIES

The subtraction method can be used to screen cDNA libraries from diverse sources. For example, the same cell types from different species can be compared by gene

15 transcript analysis to screen for specific differences, such as in detoxification enzyme systems. Such testing aids in the selection and validation of an animal model for the commercial purpose of drug screening or toxicological testing of drugs intended for human or animal use. When

20 the comparison between animals of different species is shown in columns for each species, we refer to this as an interspecies comparison, or zoo blot.

Embodiments of this invention may employ databases such as those written using the FoxBASE programming

25 language commercially available from Microsoft Corporation. Other embodiments of the invention employ other databases, such as a random peptide database, a polymer database, a synthetic oligomer database, or a oligonucleotide database of the type described in U.S. Patent 5,270,170, issued

30 December 14, 1993 to Cull, et al., PCT International Application Publication No. WO 9322684, published November 11, 1993, PCT International Application Publication No. WO 9306121, published April 1, 1993, or PCT International Application Publication No. WO 9119818, published December 26, 1991. These four references (whose text is incorporated herein by reference) include teaching which

PCT/US95/01160 WO 95/20681

5

may be applied in implementing such other embodiments of the present invention.

All references referred to in the preceding text are hereby expressly incorporated by reference herein.

Various modifications and variations of the described method and system of the invention will be apparent to those skilled in the art without departing from the scope and spirit of the invention. Although the invention has been described in connection with specific preferred 10 embodiments, it should be understood that the invention as claimed should not be unduly limited to such specific embodiments.

	•		
		TABLE 1	
Designations (D)	Distribution (F)	Localization (Z)	Function (R)
1	Non-specific	Nuclear	T = Translation
n			L = Protein processing
H = Homologous	P = Cell/clasue specialc	Z CYCOPICSINIO	R = Ribosomal protein
n	U = Unknown		O = Oncodene
·N = No match	,	THE COLLEGE BORD	
D = Noncoding gene			v = viral element
U = Nonreadable		M = MICOCHOINTIAL	v - Kinase/nhosphatase
R = Repetitive DNA		S = Secreted	
	Species	U = Unknown	
V = Vector only	(8)	X = Other	I = Binding processes
M = Mitochondrial DNA			n
4	H = Human		B = Surface molecule/receptor
7 I w Match Incyte clone	A = Ape		C = Ca binding protein
×	P = Pig		
	D = Dog		
T. Street	V - Bovine	Status	E = Enzyme
	R = Rabbit	(H)	F = Perroprotein
(4)	1		P = Protease/inhibitor
2001 - 11	MOUSE MOUSE	0 = No current interest	Z = Oxidative phosphorylation
		1 - Do orimary analysis	Q = Sugar metabolism
= HA		o = primary analysis done	M = Amino acid metabolism
T = THP-1	C = Culcken	Occopy of the state of the stat	He Note of the Challen I N
H = HUVEC	F = Amphibian	3 = Full length sequence	
S = Spleen	<pre>I = Invertebrate</pre>	4 - Secondary analysis	8
L = Lung	Z = Protozoan	5 = Tissue northern	K = Structural
Y = T & B cell	G = Fungi	6 = Obtain full length	X = Other
A - Adenoid			U = unknown

TABLE 2

Clone numbers 15000 through 20000 Libraries: HUVEC Arranged by ABUNDANCE Total clones analyzed: 5000

319 genes, for a total of 1713 Clones

	number	N c	:	entry	S	descriptor
1	15365	67		HSRPL41		Riboptn L41
2	15004	65		NCY015004		INCYTE 015004
3	15638	63		NCY015638		INCYTE 015638
4	15390	50		NCY015390		INCYTE 015390
5	15193	47		HSFIB1		Fibronectin
6	15220	47		RRRPL9	R	Riboptn L9
7	15280	47		NCY015280		INCYTE 015280
8	15583	33		M62060		EST HHCHO9 (IGR)
9	15662	31		HSACTCGR		Actin, gamma .
10	15026	29		NCY015026		INCYTE 015026
11	15279	24		HSEF1AR		Elf 1-alpha
12	15027	23		NCY015027		INCYTE 015027
13	15033	20		NCY015033		INCYTE 015033
14	15198	20		NCY015198		INCYTE 015198
15	15809	20		HSCOLL1		Collagenase
16	15221	19		NCY015221		INCYTE 015221
17	15263	19		NCY015263		INCYTE 015263
18	15290	19		NCY015290		INCYTE 015290
19	15350	18		NCY015350		INCYTE 015350
20	15030	17		NCY015030		INCYTE 015030
21	15234	17		NCY015234		INCYTE 015234
22	15459	16		NCY015459		INCYTE 015459
23	15353	15		NCY015353		INCYTE 015353
24	15378	15		S76965		Ptn kinase inhib
25	15255	14		HUMTHYB4		Thymosin beta-4
26 .	15401	14		HSLIPCR		Lipocortin I
27	15425	14		HSPOLYAB		Poly-A bp
28	18212	14		HUMTHYMA		Thymosin, alpha
29	18216	14		HSMRP1		Motility relat ptn; MRP-1;CD-9
30	15189	13		HS18D		Interferon induc ptn 1-8D
31	15031	12		HUMFKBP	_	FK506 bp
32	15306	12		HSH2AZ		Histone H2A
.33	15621	12		HUMLEC		Lectin, B-galbp, 14kDa
34	15789	11		NCY015789		INCYTE 015789
35	16578	11		HSRPS11		Riboptn S11
36	16632 ·	11		M61984		EST HHCA13 (IGR)
37	18314	11		NCY018314		INCYTE 018314
38	15367	10		NCY015367		INCYTE 015367
39	15415	10		HSIFNIN1		interferon induc mRNA
40.	15633	10		HSLDHAR		Lactate dehydrogenase
41	15813	10		CHKNMHCB .		C Myosin heavy chain B
42	18210	10		NCY018210		INCYTE 018210
43	18233	10		HSRPII140		RNA polymerase II
44	18996	10		NCY018996		INCYTE 018996
45	15088	9.		HUMFERL		Ferritin, light chain
46	15714	9		NCY015714		INCYTE 015714
47	15720	9		NCY015720		INCYTE 015720
48	15863	9		NCY015863		INCYTE 015863
49	16121	9		HSET		Endothelin
50	18252	9		NCY018252		INCYTE 018252
51	15351	8		HUMALBP		Lipid bp, adipocyte
52	15370	8		NCY015370		INCYTE 015370

TABLE 2 Con't

	number	N	C	entry	s	descriptor
53	15670	. 8		BTCIASHI	v	NADH-ubig oxidoreductase
54	15795	8		NCY015795		INCYTE 015795
55	16245	8		NCY016245		INCYTE 016245
56	18262	8		NCY018262		INCYTE 018262
57	18321	8		HSRPL17		Riboptn L17
58	15126	7		XLRPL1BRF		Riboptn L1
59	15133	7		HSAC07		Actin, beta
60	15245	7		NCY015245		INCYTE 015245
61	15288	7		NCY015288		INCYTE 015288
62	15294	7		HSGAPDR		G-3-PD
63	15442	7		HUMLAMB		Laminin receptor, 54kDa
64	15485	7		HSNGMRNA		Uracil DNA glycosylase
65	16646	7		NCY016646		INCYTE 016646
66	18003	7		HUMPAIA		Plsmnogen activ gene
67	15032	6		HUMUB		Ubiquitin
68	15267	6		HSRPS8		Riboptn S8
69	15295	6		NCY015295		INCYTE 015295
70	15458	6		RNRPS10R	R	Riboptn S10
71	15832	6		RSGALEM	R	UDP-galactose epimerase
72	15928	6		HUMAPOJ		Apolipoptn J
73	16598	6		HUMTBBM40		Tubulin, beta
74	18218	6		NCY018218		INCYTE 018218
75	18499	6		HSP27		Hydrophobic ptn p27
76	18963	6		NCY018963		INCYTE 018963
77	18997	6		NCY018997		INCYTE 018997
78	15432	5		HSAGALAR		Galactosidase A, alpha
79	15475	5		NCY015475		INCYTE 015475
80	15721	5		NCY015721		INCYTE 015721
81	15865	5		NCY015865		INCYTE 015865
82	16270	5		NCY016270		INCYTE 016270
83	16886	5		NCY016886		INCYTE 016886
84	18500	5		NCY018500		INCYTE 018500
85	18503	5		NCY018503	_	INCYTE 018503
86	19672	5		RRRPL34	R	Riboptn L34
87	15086	4		XLRPL1AR	F	Riboptn Lla
88	15113	4		HUMIFNWRS		tRNA synthetase, trp
89	15242	4		NCY015242		INCYTE 015242
90	15249	4		NCY015249		INCYTE 015249
91	15377	4		NCY015377		INCYTE 015377
92	15407	4		NCY015407		INCYTE 015407
93	15473	4 .		NCY015473		INCYTE 015473
94	15588	4		HSRPS12		Riboptn S12
95	15684	. 4		HSEF1G	•	Elf 1-gamma
96	15782	4		NCY015782		INCYTE 015782
97	15916	4		HSRPS18		Riboptn S18
98	15930	4		NCY015930		INCYTE 015930
99	16108	4		NCY016108		INCYTE 016108
100	16133	4		NCY016133		INCYTE 016133

NORMAL MONONCYTE VS. ACTIVATED MACROPHAGE

Top 15 Most Abundant Genes

NORMAL

ACTIVATED

NGF-related B cell activation molecule Macrophage inflammatory protein-l Adenylate cyclase (yeast homolog) Protease Nexin-I, glial-derived **Tumor Necrosis Factor-alpha** Rantes T-cell specific protein Cu/Zn superoxide dismutase Lymphocyte activation gene Osteopontin; nephropontin Elongation factor-l alpha Poly A binding protein NCYTE clone 01 1050 Interleukin-I beta nterleukin-8 Beta actin ranslationally controlled tumor ptn Ribosomal protein S20 homolog Signal recognition particle SRP9 Histone H2A.Z Ribosomal protein 58 homolog **Ribosomal phosphoprotein Elongation factor-I alpha** Ribosomal protein Ke-3 **Ribosomal protein S25** oly-A binding protein Ribosomal protein L7 ransferrin receptor erritin H chain **Aucleoplasmin** Beta-Globin

TABLE 3

PCT/US95/01160 WO 95/20681

TABLE 4

Libraries: THP-1 Subtracting: HMC Sorted by ABUNDANCE Total clones analyzed: 7375

1057 genes, for a total of 2151 clones

100, 90,					
number	entry s	descriptor	bgfreq	rfend	ratio
10022	HUMIL1	IL 1-beta	0	131	262.00
	HSMDNCF	IL-8	Ö	119	238.00
10036		Lymphocyte activ gene	ŏ	71	142.00
10089	HSLAG1CDN		Ö	23	46.000
10060	HUMTCSM	RANTES			
10003	HUMMIP1A	MIP-1	3	121	40.333
10689	HSOP	Osteopontin	0	20	40.000
11050	NCY011050	INCYTE 011050	0	17	34.000
10937	HSTNFR	TNF-alpha	0	17	34.000
10176	HSSOD	Superoxide dismutase	0	14	28.000
10886	HSCDW40	B-cell activ, NGF-relat	0	10	20.000
10186	HUMAPR	Early resp PMA-induc	0	9	18.000
10967	HUMGDN	PN-1, glial-deriv	0	9	18.000
11353	NCY011353	INCYTE 011353	0	8	16.000
10298	NCY010298	INCYTE 010298	0	7	14.000
10215	HUM4COLA	Collagenase, type IV	0	6	12.000
10276	NCY010276	INCYTE 010276	0	6	12.000
10488	NCY010488	INCYTE 010488	0	6	12.000
11138	NCY011138	INCYTE 011138	0	6	12.000
10037	HUMCAPPRO	Adenylate cyclase	1	10	10.000
10840	HUMADCY	Adenylate cyclase	0	5 .	10.000
10672	HSCD44E	Cell adhesion glptn	0	5	10.000
12837	HUMCYCLOX	Cyclooxygenase-2	0	5	10.000
10001	NCY010001	INCYTE 010001	0	5	10.000
10005	NCY010005	INCYTE 010005	0	5	10.000
10294	NCY010294	INCYTE 010294	0	5	10.000
10297	NCY010297	INCYTE 010297	0	5	10.000
10403	NCY010403	INCYTE 010403	0	5	10.000
10699	NCY010699	INCYTE 010699	0	• 5	10.000
10966	NCY010966	INCYTE 010966	0	5	10.000
12092	NCY012092	INCYTE 012092	0	5	10.000
12549	HSRHOB	Oncogene rho	0	5	10.000
10691	HUMARF1BA	ADP-ribosylation fctr	0	4	8.000
12106	HSADSS	Adenylosuccinate synthetase	0	4	8.000
10194	HSCATHL	Cathepsin L	0	- 4	8.000
10479	CLMCYCA I	Cyclin A	0	4	8.000
10031	NCY010031	INCYTE 010031	0	4	8.000
10203	NCY010203	INCYTE 010203	0	4	8.000
10288	NCY010288	INCYTE 010288	0	4	8.000
10372	NCY010372	INCYTE 010372	Ō	4	8.000
10471	NCY010471	INCYTE 010471	0	4	8.000
10484	NCY010484	INCYTE 010484	Ō	4	8.000
10859	NCY010859	INCYTE 010859	Ō	4	8.000
10890	NCY010890	INCYTE 010890	Ō	4	8.000
11511	NCY011511	INCYTE 011511	Ö	4	8.000
11868	NCY011868	INCYTE 011868	Ö	4	8.000
12820	NCY012820	INCYTE 012820	ŏ	4	8.000
10133	HSI1RAP	IL-1 antagonist	ŏ	4	8.000
	HUMP2A	Phosphatase, regul 2A	ŏ	4	8.000
11063	HUMB94	TNF-induc response	Ŏ	4	8.000
11140	HSHB15RNA	HB15 gene; new Ig	Ö	3	
		INCYTE 001713	ŏ	3	6.000
10788	NCY001713 NCY010033	INCYTE 010033	0	3 3	6.000
10033			0	ა 2	6.000
10035	NCY010035	INCYTE 010035		3	6.000
10084	NCY010084	INCYTE 010084	0	3	6.000
10236	NCY010236	INCYTE 010236	0	3	6.000
10383	NCY010383	INCYTE 010383	0	3	6.000

PCT/US95/01160

TABLE 4 Con't

number	entry	8	descriptor	bgfreq	rfend	ratio
10450	NCY010450		INCYTE 010450	0	3	6.000
10470	NCY010470		INCYTE 010470	0	3	6.000
10504	NCY010504		INCYTE 010504	0	3	6.000
10507	NCY010507		INCYTE 010507	0	3	6.000
10598	NCY010598		INCYTE 010598	0	3 3	6.000
10779	NCY010779		INCYTE 010779	0	3	6.000
10909	NCY010909		INCYTE 010909	0	3	6.000
10976	NCY010976		INCYTE 010976	0	3	6.000
10985	NCY010985		INCYTE 010985	0	3	6.000
11052	NCY011052		INCYTE 011052	0	3 3 3	6.000
11068	NCY011068		INCYTE 011068	0	3	6.000
11134	NCY011134		INCYTE 011134	0	3 3	6.000
11136	NCY011136		INCYTE 011136	. 0	3	6.000
11191	NCY011191		INCYTE 011191	0	3 3	6.000
11219	NCY011219		INCYTE 011219	0	3	6.000
11386	NCY011386		INCYTE 011386	0	3	6.000
11403	NCY011403		INCYTE 011403	0	3	6.000
11460	NCY011460		INCYTE 011460	0	3	6.000
11618	NCY011618		INCYTE 011618	0	3	6.000
11686	NCY011686		INCYTE 011686	0	3 3	6.000
12021	NCY012021		INCYTE 012021	0	3	6.000
12025	NCY012025		INCYTE 012025	0	3	6.000
12320	NCY012320		INCYTE 012320	0	3	6.000
12330	NCY012330		INCYTE 012330	0	3	6.000
12853	NCY012853		INCYTE 012853	0	3	6.000
14386	NCY014386		INCYTE 014386	0	3 3	6.000
14391	NCY014391		INCYTE 014391	0	3	6.000

TABLE 5

```
    Master menu for SUBTRACTION output

   BET TALK OFF
   SET SAFETY OFF
   SET EXACT ON
   SET TYPEAHEAD TO 0
   CLEAR
   SET DEVICE TO SCREEN
   USE . * SmartGuy: FoxBASE+/Mac:fox files:Clones.dbf*
   GO TOP
   STORS NUMBER TO INITIATE
   GO BOTTOM
   STORE NUMBER TO TERMINATE
                                                                         TO Target1
   STORE '
   STORE '
                                                                   ' TO Target2
                                                                    ' TO Target3
   STORE '
                                                                   ' TO Object1
   STORE.'
   STORE '
                                                                         TO Object2
   STORE '
                                                                   1 TO Object3
   STORE O TO ANAL
   STORE O TO EMATCH
   STORE O TO HMATCH
   STORE D TO OMATCH
   STORE 0 TO IMATCH
   STORE 0 TO PTY
   STORE 1 TO BAIL
   DO WHILE .T.
    * Program .: Subtraction 2.fmt
  * Date...: 10/11/94
* Version,: FooEASE+/Mac, revision 1.10
    Notes ... : Pormat file Subtraction 2
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FORT "Geneva",9 COLOR 0,0,0,0 PIXELS 75,120 TO 178,241 STYLE 3871 COLOR 0,0,-1,24610,-1,8947
6 FIXELS 27,134 SAY "Subtraction Memu" STYLE 65536 FORT "Geneva",274 COLOR 0,0,-1,-1,-1,-1 FIXELS 117,126 GET EMATCH STYLE 65536 FORT "Chicago",12 FICTURE "0"C Homologous" SIZE 15,62 CO FIXELS 135,126 GET HMATCH STYLE 65536 FORT "Chicago",12 FICTURE "0"C Homologous" SIZE 15,1 PIXELS 153,126 GET CMATCH STYLE 65536 FORT "Chicago",12 FICTURE "0"C Homologous" SIZE 15,8 FIXELS 90,152 SAY "Matches: STYLE 65536 FORT "Chicago",12 FICTURE "0"C Other spc" SIZE 15,8 FIXELS 252,137 GET initiate STYLE 65536 FORT "Chicago",12 FICTURE "0"C Incyte SIZE 15,65 CO FIXELS 252,236 GET terminate STYLE 0 FORT "Geneva",12 COLOR 0,0,-1,-1,-1,-1 FIXELS 252,235 SAY "Include clones" STYLE 65536 FORT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">" STYLE 65536 FORT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">" STYLE 65536 FORT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">" STYLE 65536 FORT "Geneva",14 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">" STYLE 65536 FORT "Geneva",14 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 252,215 SAY ">" STYLE 65536 FORT "Geneva",14 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 90,9 TO 181,109 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1 FIXELS 90,288 TO 181,397 STYLE 65536 FORT "Geneva",270 COLOR 0,0,-1,-1,-1,-1,-1 FIXELS 45,135 CET ANAL STYLE 65536 FORT "Chicago",12 FICTURE "8"R OVERALL; Function" SIZE 4
 ### PIXELS 81,296 SAY "Background:" STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1

### PIXELS 45,135 GET ANAL STYLE 65536 FONT "Chicago",12 PICTURE "%" Overall; Function" SIZE 4

### PIXELS 81,26 SAY "Target:" STYLE 65536 FONT "Geneva",270 COLOR 0,0,-1,-1,-1,-1

### PIXELS 108,20 GET target1 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 135,20 GET target2 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 108,299 GET target3 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 108,299 GET object1 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 135,299 GET object2 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 162,299 GET object3 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 162,299 GET object3 STYLE 0 FONT "Geneva",9 SIZE 12,79 COLOR 0,0,-1,-1,-1,-1

### PIXELS 276,324 GET Bail STYLE 65536 FONT "Chicago",12 PICTURE "% ### Run;Bail out" SIZE 4112
    * EOF: Subtraction.2.fmt
   READ
          IF Bail=2
           CLEAR
           CLOSE DATABASES
           USE 'SmartGuy: FoxBASE+/Mac:fox files:clones.dbf'
          SET SAFETY ON
           SCREEN 1 OFF
           RETURN
```

```
ENDIP
STORE VAL(SYS(2)) TO STARTIME
STORE UPPER(Target1) TO Target1
STORE UPPER(Target2) TO Target2
STORE UPPER(Targ t3) TO Target3
STORE UPPER(Object1) TO Object1
STORE UPPER(Object2) TO Object3
STORE UPPER(Object3) TO Object3
clear
SET TALK ON
GAP = TERMINATE-INITIATE+1
GO INITIATE
COPY NEXT GAP FIELDS NUMBER, Library, D. F. Z. R. ENTRY, S. DESCRIPTOR, START, RFEND, I TO TEMPNUM
use tempnum
COUNT TO TOT
COPY TO TEMPRED FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='I'
USE TEMPRED
IP Bratch=0 , AND, Pratch=0 , AND. Onatch=0 .AND. IMPICH=0
COPY TO TEMPDESIG
else
COPY STRUCTURE TO TEMPDESIGUSE TEMPDESIG
   IF Ematchel
   APPEND FROM TEMPNUM FOR D='B'
   ENDIF
   IF Hmatch=1
   APPEND PROM TEMPNUM FOR D='H'
   ENDIR
   IF Ometch=1
  APPEND FROM TEMPNUM FOR D='0'
  ENDIF
  IF Imatch=1
  APPEND FROM TEMPNUM FOR D='I'.OR.D='X'
*., OR, D='N'
 , EXDIF
EXILIP
COUNT TO STARTOT
COPY STRUCTURE TO TEMPLIE
.USE TEMPLIE .
  APPEND FROM TEMPDESIG FOR library=UPPER(target1)
  IP target2co'
  APPEND FROM TEMPDESIG FOR library=UPPER(target2)
  endif
  IF target3<>
  APPEND FROM TEMPDESIG FOR library=UPPER(target3)
  ENDIF
COUNT TO ANALITOT
USE TEMPDESIG
COPY STRUCTURE TO TEMPSUB
USE TEMPSUB
  APPEND FROM TEMPDESIG FOR library=UPPER(Object1)
  IF target20'
  APPEND FROM TEMPDESIG FOR library=UFFER(Object2)
  ENDIP
  IF target30'
  APPEND FROM TEMPORSIG FOR library=UPPER(Object3)
  ENDIP
COUNT TO SUBTRACTOT
BRT TALK OFF
* COMPRESSION SUBROUTINE A ? 'COMPRESSING QUERY LIBRARY'
USE TEMPLIE
```

```
SORT ON ENTRY, NUMBER TO LIBSORT
 UBE LIBSORT
 COUNT TO IDGENE
 REPLACE ALL REEND WITH 1
 MARKI - 1
 SW2=0
 DO WHILE SW2=0 ROLL
   IP MARKI >= IDGENE
PACK
   COUNT TO AUNIQUE
   SW2=1
   LOOP
   ENDIP
GO MARKL
DUF = 1
STORE ENTRY TO TESTA
STORE D TO DESIGA
 BW = 0
 DO WHILE SW=0 TEST
 EKIP
 STORE ENTRY TO TESTS
STORE D TO DESIGE
   IF TESTA = TESTB.AND.DESIGA=DESIGB
DELETE
   DUP = DUP+1
   LOOP
   EVDIF
GO MARKI
REPLACE RPEND WITH DUP
MARK1 = MARK1+DUP "
 SW=1
LOOP
ENDDO. TEST
LOOP
ENDDO ROLL
SORT ON REEND/D, NUMBER TO TEMPTARSORT.
USE TEMPTARSORT
*REPLACE ALL START WITH RPEND/IDCENE*10000
COUNT TO TEMPTARCO
• COMPRESSION SUBROUTINE B
? 'COMPRESSING TARGET LIBRARY'
USE TEMPSUB
SORT ON ENTRY, NUMBER TO SUBSORT
USE SUBSORT
COUNT TO SUBGENE
REPLACE ALL REPLO WITH 1
MARK! = 1
SW2=0
DO WHILE SW2=0 ROLL
   IF MARK1 >= SUEGENE
   PACK .
  COUNT TO BUNIQUE
  SW2=1
  LOOP
  ENDIF
GO MARKI .
DUP = 1
STORE ENTRY TO TESTA
STORE D TO DESIGA
SW = 0
DO WHILE SW=0 TEST
BKIP
STORE EXTRY TO TESTE
STORE D TO DESIGE
  IF TESTA = TESTB.AND.DESIGA=DESIGB
```

```
DELETE
DUP = DUP+1
     LOOP
ENDIP
  GO MARKI
  REPLACE REEND WITH DUP
MARK! = MARK!+DUP
  5W-1
  LOOP
 ENDO TEST
LOOF:
ENDO ROLL
SORT ON REEND/D, NUMBER TO TEMPSUBSORT
 ·USE TEMPSUBSORT
 *REPLACE ALL START WITH RPEND/LIGENE+10000
COUNT TO TEMPSUBCO
  *FUSION ROUTINE
 ? 'SUBTRACTING LIBRARIES'
USE SUBTRACTION
COPY STRUCTURE TO CRENCHER
 SELECT 2
USB TEMPSUBSORT
 SELECT 1
 USE CRINCHER
APPEND FROM TEMPTARSORT
COUNT TO BAILOUT
 MARK = 0
DO WHILE .T.
SELECT 1
MARK = MARK+1
IF MARK>BALLOUT
    EXIT
    ENDIP
 GO MARK
 STORE ENTRY TO SCANNER
SELECT 2
LOCATE FOR ENTRY-SCANNER
IP FOUND()
STORE REEND TO BIT1
STORE REEND TO BIT3
 RISE
STORE 1/2 TO BITI
STORE 0 TO BIT2
ENDIF
SELECT 1
REPLACE BGFRBQ WITH BIT2
REPLACE ACTUAL WITH BIT1
LOOP
ENDOO
SELECT 1
REFLACE ALL RATIO WITH RFEND/ACTUAL
7 'DOING FINAL SORT BY RATIO'
SORT ON RATIO/D, EGFREQ/D, DESCRIPTOR TO FINAL
USE FINAL
set talk off
DO CASE
CASE PTF=0
SET DEVICE TO PRINT
SET PRINT ON
EJECT
CASE PTF=1
SET ALTERNATE TO "Adenoid Patent Figures: Subtraction.txt"
```

ENDIP

```
SET ALTERNATE ON
ENICASE
STORE VAL(SYS(2)) TO FINTIME
IF FINTIME-STARTIME
STORE PINTING+86400 TO PINTING
ENDIP
STORE FINTING - STARTING TO COMPSEC
STORE COMPSEC/60 TO COMPMIN
*****
SET MARGIN TO 10
81,1 SAY "Library Subtraction Analysis" STYLE 65536 FONT "Geneva",274 COLOR 0,0,0,-1,-1,-1
? date()
?? TIME()
?? Clone numbers
?? SIR(INITIATE, 5, 0)
?? SIR(TEMINATE, 6, 0)
?! Libraries:
? Target1
IP Target2<>'
77 Target2
EVOLF
IP Target3<>'
?? ',
?? Target3
ENDI?
? Subtracting:
 ? Object1
 IF-Object2<>
?? Object2
ENDIF
 IF Object3<>'
Properties

?? 'Djectl

PMOT:

? 'Designations:'

IF Smatch=0 .AND. Hmatch=0 .AND. Omatch=0 .AND. IMATCH=0

?? 'All'
ENDIF ..
IF Ematch=1
?? 'Exact,'
ENDIF
IF Hmatch=1
?? 'Human,'
ENDIF .
'IF Ometch=1
?? 'Other sp.'
ENDIP
IF Datch=1
?? 'INCYTE'
ENDIF
IF ANAL=1
? 'Sorted by ABUNDANCE'
ENDLY.
IF ANAL=2
 ? 'Arranged by FUNCTION'
```

```
? 'Total clones represented: '
 ?? STR (TOT, 5, 0)
 ? 'Total clones analyzed: '
 ?? STR (STARTOT, 5, 0)
? 'Total computation time:

?? STR (COMPMIN, 5, 2)

?? ' minutes'
 ?''d = designation f = distribution s = location <math>r = function
                                                                                               s = species
                                                                                                                  i = inte
 *****
 SCREEN 1 TYPE 0 HEADING Soreen 1º AT 40,2 SIZE 286,492 PIXEL6 FONT Geneva 9 COLOR 0,0,0,
 DO CASE
 CASE ANAL-1
 ?? STR(AUNIQUE,4,0)
?? 'genes, for a total of '
?? STR (ANALMOT, A, 0)
?? 'clones'
 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0.0.0.
 list OFF fields number, D. F. Z. R. ENTRY, S. DESCRIPTOR, EGFREQ, RFEND, RATIO, I
 SET PRINT OFF
 CLOSE DATABASES
USE 'SmartGuy: FoxBASE+/Mac: fox files: clones.dbf'
CASE. ANAL=2
· * arrange/function
SET PRINT ON
 SET HEADING ON
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",268 COLOR 0
                                          BINDING PROTEINS
SCREEN 1 TYPE 0 READING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR 0
7 'Surfada molecules and receptors:'
SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Geneva',7 COLOR 0,0,0,
list CFF fields number,D,F;Z,R,EVTRY,S,DESCRIPTOR,BGFREQ,RFEND,RATTO,I FOR R='B'
SCREWN 1 TYPE 0 HEADING "Screen 1" AT 40.2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
9 'Calcium-binding proteins:'
SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Geneva',7 COLCR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATTO, I FOR R='C'
 SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica".265 COLOR O
? 'Ligands end effectors:'
SCREEN 1 TYPE 0 HEALING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
list OFF fields number,D,F,Z,R,ENTRY,S,DESCRIPTOR,BGFREQ,RFEND,RATIO,I FOR R='5'
 SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR O
 ? 'Other binding proteins:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFRSQ, RFEND, RATIO, I FOR R='I'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",268 COLOR 0 CNCCGENES"
 SCREEN 1 TYPE O HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Helvetica',265 COLOR O
 ? 'General oncogenes:'
7 'General oncogenes:'.
SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Geneva',7 COLOR 0,0,0,
list OFF fields number, D, F, Z, R, ENTRY, 8, DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR R='0'
 SCREEN 1 TYPE O HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FONT 'Helvetica',265 COLOR O
 7 'GTP-binding proteins!'
SCREEN 1 TYPE 0 HEADING 'Screen 1' AT 40,2 SIZE 286,492 PIXELS FORT 'Geneva',7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFRED, RFEND, RATIO, I FOR R='G'
```

```
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR 0
 ? 'Viral elements!'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FORT -CERES , 7 COLOR 3,0,0,
 list OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR REVY
 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0
 ? 'Kinases and Phosphatases!' .
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0, list OFF fields number, D.F. Z. R. ENTRY, S. DESCRIPTOR, BGFRED, RFEND, RATIO, I FOR R='Y'
 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PINELS FORT "Helvetica", 265 COLOR 0
 ? 'Timor-related antigens:'
SCREEN 1 TYPE D HEADING "Screen 1" AT 40.2 SIZE 286;492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number,D,F,Z,R,ENTRY,S,DESCRIPTOR,BGFREQ,RFEND,RATIO,I FOR R='A'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",268 COLOR 0
                      PROTEIN SYNTHETIC MACHINERY PROTEINS!
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZ3 286,492 PIXELS FONT "Helvetica".265 COLOR 0 7 "Transcription and Nucleic Acid-binding proteins:"
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number,D,F,Z,R,ENTRY.S,DESCRIPTOR,BGFREQ,RFEND,RATIO,I FOR R='D'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FONT "Helvetica",265 COLOR 0
 ? 'Translation:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, EGFREQ, RFEND, RATIO, I FOR R="T"
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FORT "Helvetica",265 COLOR 0
? 'Ribosomal proteins:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number.D;F,Z,R,ENTRY,S,LESCRIPTOR,BGFREQ,RFEND,RATIO,I FOR R='R'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR 0
? 'Protein processing:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
list OFF fields number, D. F. Z. R. ENTRY, S. DESCRIPTOR, EGFREQ, RFEND, RATTO, I FOR Re'L'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Helvetica",268 COLOR 0
                                           ENZYMES!
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR 0
? 'Ferroproteins:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FONT "Geneva",7 COLOR 0.0,0, list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RPEND, RATIO, I FOR R='F'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 285,492 FIXELS FORT "Helvetica",265 COLOR 0
? 'Proteases and inhibitors:'
SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
list OFF fields number, D, P, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R= 'p'
SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FORT "Helvetica",265 COLOR O
7 'Oxidative phosphorylation:'..

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,0 list OFF fields number, D, F, Z, R, ENTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='Z'
SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Helvetica",265 COLOR O
7 'Sugar metabolism:'
SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
list OFF fields number, D, F, Z, R, EMTRY, S, DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='Q'
```

SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 285,492 PIXELS FORT "Helvetica",265 COLOR 0

SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,

7 'Amino acid metabolism:

PCT/US95/01160

list OFF fields number, D.F. E.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='M'

SCREEN 1 TYPE 0. HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Helvetica", 265 COLOR 0

? 'Nucleic acid metabolism;'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Geneva", 7 COLOR 0,0,0,
list OFF fields number, D.F. Z.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='N'

ECREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Helvetica", 265 COLOR 0

? 'Lipid metabolism:'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Geneva", 7 COLOR 0,0,0,
list OFF fields number, D.F. Z.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='W'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Helvetica", 265 COLOR 0

? 'Other scrymes:'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Geneva", 7 COLOR 0,0,0,
list OFF fields number, D.F. Z.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='E'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Geneva", 7 COLOR 0,0,0,
list OFF fields number, D.F. Z.R. ENTRY, S. DESCRIPTOR, BGFREQ, RFEND, RATIO, I FOR R='E'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FCNT "Helvetica", 268 COLOR 0

**MISCELLANEOUS CATEGORIES'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 ? 'Stress response:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number,D,F,Z,R,ENTRY,S,DESCRIPTOR,BGFREQ,RFEND,RATIO,I FOR R='H'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR'O ? 'Structural:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number,D,F,Z,R,ENTRY.8,DESCRIPTOR,EGFREQ,RFEND,RATIO,I FOR R='K'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 ? 'Other clones: 'SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0 list OFF fields number,D,P,Z,R,ENTRY,S,DESCRIPTOR,BGFREQ,RPEND,RATIO,I FOR R='X'

SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Helvetica",265 COLOR 0 ? 'Clones of unknown function:'
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0, list OFF fields number,D,F,Z,R,ENIRY,S,DESCRIPTOR,EGFREQ,RFEND,RATIO,I FOR R='U'
ENDCASE

DO "Test print.pry"
SET FRINT OFF
SET DEVICE TO SCREEN
CLOSE DATABASES
ERASE TEMPLES. DEF
ERASE TEMPLES. DEF
ERASE TEMPLESIG. DEF
SET MARGIN TO 0
CLEAR
LOOP
ENDOO

5 5

```
"Northern (single), version 11-25-94
      close databases
     SET TALK OFF
SET PRINT OFF
     SET EXACT OFF
     CLEAR .
     STORE !
                                             ' TO Eobject
     STORE
                                                                                            ' TO Dobject
    STORE 0 TO Numb store 0 TO Zog
STORE 1 TO Bail
    DO WHILE .T.
    * Program. : Northern (single) .fmt
     * Data...: 8/ 8/94
    * Version .: . FoxBASE+/Mac, revision 1.10
    * Notes.:..: Format file Northern (single)
 SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",12 COLOR 0,0,0 PIXELS 15,81 TO 46,397 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1
6 PIXELS 89,79 TO 192,422 STYLE 28447 COLOR 0,0,0,-25600,-1,-1
6 PIXELS 89,79 TO 192,422 STYLE 28447 COLOR 0,0,0,-25600,-1,-1
6 PIXELS 115,98 SAY "Entry %: STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
7 PIXELS 15,173 GET Expiect STYLE 0 FONT "Geneva",12 SIZE 15,142 COLOR 0,0,0,-1,-1,-1
8 PIXELS 145,89 SAY "Description" STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
8 PIXELS 145,173 GET Dobject STYLE 0 FONT "Geneva",12 SIZE 15,241 COLOR 0,0,0,-1,-1,-1
8 PIXELS 35,89 SAY "Single Morthern search screen" STYLE 65536 FONT "Geneva",274 COLOR 0,0,-
8 PIXELS 220,162 GET Bail STYLE 65536 FONT "Chicago",12 PICTURE "@*R Continue;Bail cut" SIZE
8 PIXELS 175,98 SAY "Clome %:" STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
8 PIXELS 175,173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,0,-1,-1,-1
8 PIXELS 175,173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,0,-1,-1,-1
8 PIXELS 80,152 SAY "Enter any CNE of the following:" STYLE 65536 FONT "Geneva",12 COLOR -1,
   * DOF: Northern (single) , fint
  READ
  IF Bail-2
  CLEAR
  screen 1 off
 RETURN
 EXDIF
 USE SmartGuy: FoxBASE+/Mac:Fox files:Lookup.dbf*
  SET TALK ON
 IF Pobject>'
STORE UPPER(Enbject) to Pobject
 SET SAFETY OFF
 SORT ON Entry TO "Lookup entry.dbf"
 SET SAFETY ON
USE 'Lookup entry dof'
LOCATE FOR LookeEobject
 IF .. NOT. FOUND()
 CLEAR
 LOOP
 PAINTE
 BROWSK
STORE Entry TO Searchval
CLOSE DATABASES
ERASE "Lookup entry dbf"
FNOTE
IF Dobjecto'.
SET EXACT OFF
SET SAPETY OFF
SORT ON descriptor TO "Lookup descriptor.dbf"
SET SAFETY On
USE 'Lookup descriptor.cbf'
LOCATE FOR UPPER(TRIM(descriptor))=UPPER(TRIM(Dobject)) .
IF .NOT.FOUND()
CLEAR
```

WO 95/20681

```
LOOP
 ENDIP
 BROWSE
 STORE Entry TO Searchval
CLOSE DATABASES
ERASE "Lookup descriptor.dbf"
 BET EXACT ON
 ENDIP .
 IF Number-0
 USE *SmartGuy: FoodASE+/Mac:Fox files:clones.dbf*
GO Munio
ERONSE
 STORE Entry TO Searchval
 ENDIF
 CLEAR
 ? 'Northern analysis for entry '
 ?? Searchval
 ? 'Enter Y to proceed'
WAIT TO OK
 CLEAR
 IF UPPER (OK) . 'Y'
 screen 1 off
RETURN
ENDIF
* COMPRESSION SUBROUTINE FOR Library.dbf
7 'Compressing the Libraries file now...'
USE 'SmartGuy:FoodBASE+/Mac:Fox files:libraries.dbf'
SET SAFETY CFF
SORT ON library TO "Compressed libraries.dbf"
* FOR entered>0
SET SAFETY ON
USE 'Compressed libraries.dbf"
DELSTE FOR entered=0
PATK
PACK
COUNT TO TOT
MARK1 = 1
 SW2=0.
DO WHILE SW2=0 ROLL
   'IF MARKI >= TOT
  PACK
5W2=1
   LOOP
   ENDIF
GO MARKI
STORE library TO TESTA
SKIP
STORE Library TO TESTE
IF TESTA = TESTE
DELETE
ENDIF
MARKI = MARKI+1
LOOP
 ENDDO ROLL
 * Northern analysis
CLEAR
7 'Doing the northern now...
SET TALK ON
USE "SmartGuy: Foo@ASE+/Mac: Fox files: clones.dbf"
 SET SAVETY OFF
COPY TO "Hits.dbf" FOR entry-searchval
 SET SAFETY ON
```

```
    MASTER ANALYSIS 3, VERSION 12-9-94

  * Master menu for analysis output
 CLOSE DATABASES
  SET TALK OFF
 SET SAFETY OFF
 CLEAR
  SET DEVICE TO SCREEN
 SET DEFAULT TO "SmartGuy: FoxBASE+/Mac: fox files:Output programs:"
 USE "SmartGuy: FoxBASE+/Mac: fox files: Clones.dbf"
 GO TOP
  STORE NUMBER TO INITIATE
 GO BOTTOM
 STORE NUMBER TO TERMINATE
STORE 0 TO ENTIRE
 STORE O TO CONDEN
 STORE O TO ANAL
 STORE 0 TO EMATCH
STORE 0 TO HMATCH
  STORE O TO OMATCH
 STORE O TO IMATCH
 STORE O TO XMATCH
 STORE O TO PRINTON
 STORE 0 TO PTF
 DO WHILE .T.
  * Program.: Master analysis.fmt
  * Date...: 12/ 9/94
 * Version.: FoxBASE+/Mac, revision 1.10
* Notes...: Format file Master analysis
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",9 COLOR 0,0,0, 6 PIXELS 39,255 TO 277,430 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1 6 PIXELS 75,120 TO 178,241 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1
PIXELS 75,120 TO 178,241 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1
PIXELS 27,98 SAY "Customized Output Menu" STYLE 65536 FCNT "Geneva",274 COLOR 0,0,-1,-1,-1
PIXELS 45,54 GET conden STYLE 65536 FONT "Chicago",12 PICTURE "@*C Condensed format" SIZE
PIXELS 54,261 GET anal STYLE 65536 FONT "Chicago",12 PICTURE "@*RV Sort/mmber; Sort/entry;
PIXELS 117,126 GET EMATCH STYLE 65536 FONT "Chicago",12 PICTURE "@*C Homologous" SIZE 15,62 CO
PIXELS 135,126 GET HMATCH STYLE 65536 FONT "Chicago",12 PICTURE "@*C Homologous" SIZE 15,1
PIXELS 153,126 GET CMATCH STYLE 65536 FONT "Chicago",12 PICTURE "@*C Other spc" SIZE 15,84
PIXELS 90,152 SAY "Matches:" STYLE 65536 FONT "Chicago",12 PICTURE "@*C Other spc" SIZE 15,84
PIXELS 90,152 GET PRINTON STYLE 65536 FONT "Chicago",12 PICTURE "@*C Include clone listing"
PIXELS 171,126 GET Imatch STYLE 65536 FONT "Chicago",12 PICTURE "@*C Incyte" SIZE 15,65 CO
PIXELS 252,146 GET initiate STYLE 0 FONT "Chicago",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1
@ PIXELS 252,146 GET initiate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1 PIXELS 270,146 GET terminate STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,-1,-1,-1,-1 PIXELS 234,134 SAY "Include clones " STYLE 65536 FONT "Geneva",12 COLOR 0,0,-1,-1,-1,-1 PIXELS 270,125 SAY "->" STYLE 65536 FONT "Geneva",14 COLOR 0,0,-1,-1,-1,-1 PIXELS 198,126 GET PIF STYLE 65536 FONT "Chicago",12 PICTURE "@*C Print to file" SIZE 15,9 PIXELS 189,0 TO 257,120 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1 PIXELS 209,8 SAY "Library selection" STYLE 65536 FONT "Geneva",266 COLOR 0,0,-1,-1,-1,-1 PIXELS 227,18 GET ENTIRE STYLE 65536 FONT "Chicago",12 PICTURE "@*RV All; Selected" SIZE 16
 * EOF: Master analysis.fmt
READ
      IF ANAL=9
      CLEAR
      CLOSE DATABASES
      ERASE TEMPMASTER.DBF
      USE "SmartGuy:FoxBASE+/Mac:fox files:clones.dbf"
      SET SAFETY ON
      SCREEN 1 OFF
      RETURN
     ENDIF
clear
 ? INITIATE
 ? TERMINATE
? . CONDEN
? ANAL
```

```
? ematch
 ? Hmatch
   Cmatch
 ? DMATCH
 SET TALK ON
   IF ENTIRE=2
 USE 'Unique libraries dbf'
   REPLACE ALL 1 WITH
   BROWSE FIELDS i, libname, library, total, entered AT 0,0
   ENDIP
 USE "SmartGuy: FoxBASE+/Mac: fox files: clones.dof"
 *COPY TO TEMPNUM FOR NUMBER>=INITIATE.AND.MCMEER<=TERMINATE
 *USE TEMPNUM
 COPY STRUCTURE TO TEMPLIB
 USE TEMPLIB
   IP ENTIRE-1
   APPEND FROM 'SmartGuy: FoxBASE+/Mac:fox files:Clones.dbf'
   ENDIF
   IF ENTIRE 2
 USE *Unique libraries.dbf*
   COPY TO SELECTED FOR UPPER (i) = 'Y'
   USE SELECTED
   STORE RECCOUNT() TO STOPIT
   MARK=1
     DO WHILE .T.
     IF MARK>STOPIT
     CLEAR .
     EXIT
     ENDIF
     USE SELECTED
     GO MARK
     STORE library TO THISONE
     ? 'COPYING '
     ?? THISONE
    USE TEMPLIE
     APPEND FROM 'SmartGuy: FoxBASE+/Mac: fox files: Clones.dbf' FOR library=THISONE
     STORE MARK+1 TO MARK
    LOOP
    ENDDO
  ENDIP
USE "SmartGuy: FoxBASE+/Mac:fox files:clones.dbf"
COUNT TO STARTOT
COPY STRUCTURE TO TEMPDESIG
USE TEMPDESIG
  IF Ematch=0 .AND. Hmatch=0 .AND. Omatch=0 .AND. IMATCH=0 APPEND FROM TEMPLIB
  ENDIP
  IF Ematch=1
  APPEND FROM TEMPLIB FOR D='E'
  ENDIF
  IF Hmatch=1
  APPEND FROM TEMPLIB FOR D='H'
  ENDIP
  IF Omatch=1
  APPEND FROM TEMPLIB FOR D='O'
  ENDIP
  IF Imatch=1
  APPEND FROM TEMPLIB FOR D='I'.OR.D='X'.OR.D='N'
  ENDIF
  IF Xmatch=1
  APPEND FROM TEMPLIB FOR D='X'
  ENDIF
COUNT TO AVALATOT
set talk off
DO CASE
```

```
CASE PTF=0
 SET DEVICE TO PRINT
 SET PRINT ON
EJECT
CASE PTF=1
 SET ALTERNATE TO "Total function sort.txt"
 *SET ALTERNATE TO "H and O function sort.txt"
*SET ALTERNATE TO "Shear Stress HOVEC 2:Abundance sort.txt"
*SET ALTERNATE TO "Shear Stress HOVEC 2:Abundance con.txt"
*SET ALTERNATE TO "Shear Stress HOVEC 2:Function sort.txt"
*SET ALTERNATE TO "Shear Stress HUVEC 2:Distribution sort.txt"
*SET ALTERNATE TO "Shear stress HUVEC 1:Clone list.txt"
 *SET ALTERNATE TO "Shear Stress HUVEC 2:Location Bort.txt"
SET ALTERNATE ON
ENDCASE
IF PRINTON=1
91,30 SAY "Database Subset Analysis" STYLE 65536 FONT "Geneva",274 COLOR 0,0,0,-1,-1,-1
ENDIF '
? date()
??'
?? TIME()
? 'Clone numbers '
?? STR(INITIATE, 6, 0)
?? 'through '
?? STR (TERMINATE, 6, 0)
? 'Libraries:
IF EVFIRE=1
? 'All libraries'
ENDIF
IF ENTIRE=2
     MARK-1
     DO WHILE .T.
     IF MARK>STOPIT
     EXIT
     ENDIF
     USE SELECTED
     GO MARK
     ?? TRIM(libname)
     STORE MARK+1 TO MARK
     LOOP
    ENDDO
ENDIF
? 'Designations: '
IF Ematch=0 .AND. Hmatch=0 .AND. Omatch=0 .AND. IMATCH=0
?? 'A11'
ENDIF
IF Ematch=1
?? 'Exact,'
ENDIP
IF Hmatch=1
?? 'Human, '
ENDIF
IF Quatch=1
?? 'Other sp. '
ENDIF
IF Imatch=1
??'INCYTE'
ENDIF
IF Xmatch=1
?? 'EST'
```

```
ENDIF
IF CONDEN-1
? 'Condensed format analysis'
ENDIF
IF ANAL=1
? 'Sorted by NUMBER'
ENDIF
IF ANAL=2
? 'Sorted by ENTRY'
ENDIF
IF ANAL=3
? 'Arranged by ABUNDANCE'
ENDIF
IF ANAL=4
? 'Sorted by INTEREST'
ENDIP
IF ANAL=5
? 'Arranged by LOCATION'
ENDIF
IF ANAL=6
? 'Arranged by DISTRIBUTION'
ENDIF
IF ANAL=7
? 'Arranged by FUNCTION'
ENDIF
? 'Total clones represented: '
?? STR(STARTOT, 6, 0)
? 'Total clones analyzed: '
?? STR(ANALITOT, 6, 0)
7 'l = library
                   d = designation f = distribution z = location r = function c = cer
******
USE TEMPDESIG
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",7 COLOR 0,0,0,
DO CASE
CASE ANAL=1
* sort/number
SET HEADING ON
IF CONDEN=1
SORT TO TEMP1 ON ENTRY, NUMBER
DO "COMPRESSION number.PRG"
ELSE
SORT TO TEMP1 ON NUMBER
USE TEMP1
list off fields number, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR
*list off fields number, L, D, F, Z, R, C, EYTRY, S, DESCRIPTOR, LENGTH, RPEND, INIT, I
CLOSE DATABASES
ERASE TEMP1.DBF
ENDIF
CASE ANAL=2
* sort/DESCRIPTOR
SET HEADING ON
*SORT TO TEMP1 ON DESCRIPTOR, ENTRY, NUMBER/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
*SORT TO TEMP1 ON ENTRY, DESCRIPTOR, NUMBER/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
SORT TO TEMP1 ON ENTRY, START/S for D='E'.OR.D='H'.OR.D='O'.OR.D='X'.OR.D='I'
IF CONDEN=1
DO "COMPRESSION entry.PRG"
elŝe
USE TEMP1
list off fields number, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, RFEND, INIT, I CLOSE DATABASES
ERASE TEMP1.DBF
ENDIF
```

```
CASE ANAL=3
 * sort by abundance
 SET HEADING ON
 SORT TO TEMP1 ON ENTRY, NUMBER for D='E'.OR.D='H'.OR.D='O'.OR.D='X',OR.D='I'
 DO "COMPRESSION abundance.PRG"
 CASE ANAL=4
 sort/interest
 SET HEADING ON
 IF CONDEN=1
 SORT TO TEMP1 ON ENTRY, NUMBER FOR I>0
 DO "COMPRESSION interest.PRG"
ELSE
 SORT ON I/D, ENTRY TO TEMP1 FOR I>1
USB TEMP1
 list off fields number, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, RFEND, INIT, I
 CLOSE DATABASES
ERASE TEMP1.DBF
ENDIP
CASE ANAL=5
 * arrange/location
SET HEADING ON
STORE 4 TO AMPLIFIER
? 'Nuclear:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I. COMMEN
IF CONDEN=1
DO "Compression location.prg"
KLSE
DO "Normal subroutine 1"
ENDIF
? 'Cytoplasmic:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
RISE
DO "Normal subroutine 1"
ENDIF
? 'Cytoskeleton:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO *Compression location.prg*
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Cell surface:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Intracellular membrane:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
PNDIP
? 'Mitochondrial:'
SORT ON ENTRY NUMBER FIELDS REED, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENOTH, INIT, I, COMMEN
IF CONDENS1
DO "Compression location.prg"
ELSE.
DO "Normal subroutine 1"
ENDIP
```

```
? 'Secreted:'
SORT ON ENTRY NUMBER FIELDS REEND, NUMBER, L.D.F.Z.R.C.ENTRY, S.DESCRIPTOR, LENGTH, INIT, I.COMMEN
IF CONDENS1
DO "Compression location.prg"
DO "Normal subroutine 1"
ENDIF
? 'Other:'
SORT ON ENTRY, NUMBER FIELDS RESD, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Unknown: '
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L. D. F. Z.R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I. COMMEN
IF CONDEN-1
DO "Compression location.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
IF CONDEN-1
SET DEVICE. TO PRINTER
SET PRINTER ON
EJECT
DO "Output heading.prg"
USE "Analysis location.dbf"
DO "Create bargraph.prg"
SET HEADING OFF
                                                           UNIQUE NEW & TOTAL'
         FUNCTIONAL CLASS
                                                   TOTAL
LIST OFF FIELDS Z, NAME, CLONES, GENES, NEW, FERCENT, GRAPH
CLOSE DATABASES
ERASE TEMP2.DBF
SET HEADING ON
*USE *SmartGuy: FoxBASE+/Mac: fox files: TEMPMASTER.dbf*
ENDIF
CASE ANAL=6
* arrange/distribution
SET HEADING ON
STORE 3 TO AMPLIFIER
? 'Cell/tissue specific distribution:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D.F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO 'Compression distrib.prg'
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Non-specific distribution:'
SORT ON ENTRY, NUMBER FIELDS RYEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression distrib.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Unknown distribution:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO "Compression distrib.prg"
RL9E
DO "Normal subroutine 1"
ENDIF
TP CONDEN=1
SET DEVICE TO PRINTER
SET PRINTER ON
```

```
EJECT
 DO "Output heading.prg"
 USE "Analysis distribution.dbf"
 DO 'Create bargraph.prg'
 SET HEADING OFF
          FUNCTIONAL CLASS
                                                    TOTAL
                                                             UNIQUE & TOTAL
 LIST OFF PIELDS P. NAME, CLONES, GENES, PERCENT, GRAPH
 CLOSE DATABASES
 ERASE TEMP2.DBF
 SET HEADING ON
 *USE *SmartGuy:FoxBASE+/Mac:fox files:TEMPWASTER.dbf*
 ENDIF
 CASE ANAL=7
 * arrange/function
 SET HEADING ON
 STORE 10 TO AMPLIFIER
                                  BINDING PROTEINS'
 ? 'Surface molecules and receptors:'
 SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
DO *Compréssion function.prg*
ELSE
DO 'Normal subroutine 1'
ENDIF
 ? 'Calcium-binding proteins:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO *Compression function.prg*
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Ligands and effectors'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIF
? 'Other binding proceins:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L.D. F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IP CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIP
*EJECT
? '
                                  ONCOGENES'
? 'General oncogenes:'
SORT ON ENTRY, NUMBER PIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'GTP-binding proteins('
SORT ON ENTRY, NUMBER FIELDS RFEND, MIMBER, L, D, P, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIP
? 'Viral elements:'
```

```
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT. I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Kinases and Phosphatases:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Tumor-related antigens:'
SORT ON ENTRY NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
*EJECT
                               PROTEIN SYNTHETIC MACHINERY PROTEINS'
 'Transcription and Nucleic Acid-binding proteins:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIP
? 'Translation:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Ribosomal proteins:'
SORT ON ENTRY, MIMBER FIELDS RFEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN-1
DO 'Compression function.prg'
ELSE
DO "Normal subroutine 1"
ENDIF
? 'Protein processing:'
SORT ON ENTRY MIMBER FIELDS RFEND, NUMBER, L, D, F, 2, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg".
DO "Normal subroutine 1"
ENDIF
*EJECT
                                  ENZYMES'
  'Ferroproteins '
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIP
? 'Proteases and inhibitors:'
SORT ON ENTRY, NUMBER FIELDS REPEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IF CONDEN=1
DO "Compression function.prg"
else
```

```
DO "Normal subroutine 1"
 ENDIP
 ? 'Oxidative phosphorylation:'
 SORT ON ENTRY, NUMBER FIELDS RPEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDENSI
DO "Compression function.prg"
 ELSE
 DO "Normal subroutine 1"
 ENDIF
 ? 'Sugar metabolism:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT; I, COMMEN
 IF CONDEN-1
 DO *Compression function.prg*
 ELSZ
 DO "Normal subroutine 1"
 ENDIF
 ? 'Amino acid metabolism:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN-1
 DO "Compression function.prg"
 E.SE
 DO "Normal subroutine 1"
 ENDIP
 ? 'Mucleic acid metabolism:'
 SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F. Z.R.C., ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
DO *Compression function.prg*
 ELSE
DO 'Normal subroutine 1"
ENDIP
 ? 'Lipid metabolism:'
SORT ON ENTRY. NUMBER FIELDS RFEND, NUMBER, L.D. F.Z.R.C. ENTRY.S. DESCRIPTOR, LENGTH, INIT, I.COMMEN
IF CONDEN=1
DO "Compression function.prg"
ELSE
DO 'Normal subroutine 1'
ENDIP
? 'Other enzymes!'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F.Z.R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMPAN
IF CONDEN=1
DO 'Compression function.prg'
ELSE
DO 'Normal subroutine 1"
ENDIF
*EJECT
                                  MISCELLANEOUS CATEGORIES'
? 'Stress'response:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
IP CONDEN-1
DO 'Compression function, prg"
ELSE
DO 'Normal subroutine 1"
ENDIP
? 'Structural:'
SORT ON ENTRY, NUMBER FIELDS REEND, NUMBER, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I, COMMEN.
IF CONDEN=1
DO 'Compression function.prg'
ELSE
DO "Normal subroutine 1"
ENDIP
? 'Other clones: '
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D.F.Z.R.C., ENTRY, S.DESCRIPTOR, LENGTH, INIT, I, COMMEN
IP CONDEN=1
DO "Compression function.prg"
ELSE
```

```
DO "Normal subroutine 1"
 ENDIP
 ? 'Clones of unknown function:'
SORT ON ENTRY, NUMBER FIELDS RFEND, NUMBER, L.D. F. Z. R.C. ENTRY, S. DESCRIPTOR, LENGTH, INIT, I, COMMEN
 IF CONDEN=1
DO "Compression function.prg"
ELSE
DO "Normal subroutine 1"
ENDIP
IF CONDEN=1
EJECT
*SET DEVICE TO PRINTER
 *SET PRINT ON
DO "Output heading.prg"
USE "Analysis function.dbf"
DO "Create bargraph.prg"
SET HEADING OFF
***
SCREEN 1 TYPE O HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",12 COLOR 0.0.0
? '
                                                                    TOTAL
                                                                               TOTAL
                                                                                                   DIST
                                                                                         NEW
? '
            FUNCTIONAL CLASS
                                                      CLONES
                                                                GENES GENES
                                                                                  FUNCTIONAL CLASS'
۱ ج٠
***
*LIST OFF FIELDS P, NAME, CLONES, GENES, NEW, PERCENT, GRAPH, COMPANY
LIST OFF FIELDS P, NAME, CLONES, GENES, NEW, PERCENT, GRAPH
CLOSE DATABASES
ERASE TEMP2.DEF
SET HEADING ON
*USE *SmartGuy: FoxBASE+/Mag: fox files: TEMPMASTER.dof*
ENDIF
CASE ANAL-8
DO "Subgroup summary 3.prg"
ENDCASE
DO "Test print.prg"
SET PRINT OFF
SET DEVICE TO SCREEN
CLOSE DATABASES
*ERASE TEMPLIB.DBP
*ERASE TEMPNUM.DBF
*ERASE TEMPDESIG.DEF
*BRASE SELECTED. DEF
CLEAR
LOOP
ENDDO
```

PCT/US95/01160

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
USE TEMP1
COUNT TO TOT
REPLACE ALL REEND WITH 1
MARK1 = 1
SW2=0
DO WHILE SW2=0 ROLL
  IF MARK1 >= TOT
  PACK
  COUNT TO UNIQUE
  COUNT TO NEWGENES FOR D='H'.OR.D='O'
  $W2=1
  LOOP
  ENDIF
GO MARKI
DUP = 1
STORE ENTRY TO TESTA
5W = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTE
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LOOP
  ENDIP
GO MARKI
REPLACE REEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
GO TOP
STORE Z TO LOC '
USE "Analysis location.dbf"
LOCATE FOR Z-LOC
REPLACE CLONES WITH TOT
REPLACE GENES WITH UNIQUE
REPLACE NEW WITH NEWGENES
USE TEMP1
SORT ON RFEND/D TO TEMP2
USE TEMP2
?? STR(UNIQUE,5,0)
?? ' genes, for a total of '
?? STR(TOT,5,0)
??'.clones'
                             v Coincidence'
list off fields number, RFEND, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE TEMPDESIG
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
USE TEMP1
COUNT TO TOT
REPLACE ALL RFEND WITH 1
MARK1 = 1
5W2=0
DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
   5W2=1
   LOOP
   ENDIP
GO MARKI
DUP = 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTB
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LOOP
 · ENDIP
GO MARKI
REPLACE RPEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON DATE TO TEMP2
USE TEMP2
?? STR(UNIQUE, 4,0)
?? ' genes, for a total of'
?? STR(TOT,4,0)
?? - clones'
                             V Coincidence'
COUNT TO P4 FOR I=4
IF P4>0
? STR(P4,3,0)
?? 'genes with priority = 4 (Secondary analysis:)'
list off fields number, RFEND, L.D.F., 2, R.C. ENTRY, S.DESCRIPTOR, LENGTH, INIT for I=4
ENDIF
COUNT TO P3 FOR I=3
IF P3>0
? STR(P3,3,0)
?? ' genes with priority = 3 (Full insert sequence:)'
list off fields number.RFEND, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR, LENGTH, INIT for I=3
ENDIP
COUNT TO P2 FOR I=2.
IF P2>0
? STR(P2,3,0)
?? ' genes with priority = 2 (Primary analysis complete:)'
list ff fields number.RFEND,L,D,F,Z,R,C,ENTRY,6,DESCRIPTOR,LENGTH,INIT for I=2
ENDIF
COUNT TO P1 FOR I=1
IF P1>0
```

? STR(P1,3,0)
?? ' genes with priority = 1 (Primary analysis needed:)'
list off fields number, RFEND, L, D, P, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT for I=1
ENDIF

SET PRINT OFF CLOSE DATABASES ERASE TEMP1.DBF ERASE TEMP2.DBF USE 'SmartGuy:FoxBASE+/Mac:fox files:clones.dbf

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 USE TEMP1
 COUNT TO TOP
REPLACE ALL RFEND WITH 1
MARK1 = 1
SW2=0
 DO WHILE SW2=0 ROLL
   IF MARKL >= TOT
   PACK
   COUNT TO UNIQUE
   5W2=1
   LOOP
   ENDIF
 GO MARK1
DUP = 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTS
IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LCOP
  EWDIP
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON NUMBER TO TEMP2
USE TEMP2
?? STR(UNIQUE,4.0)
?? 'genes, for a total of '
?? STR(TOT,5.0)
?? 'clones'
                            V Coincidence'
list off fields number.RFEND,L,D,F,Z,R,C,EMTRY,S,DESCRIPTOR,LENGTH,INIT,I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE 'SmartGuy: FoxBASE+/Mac:fox files:clones.dof'
```

PCT/US95/01160

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
USE TEMP1
COUNT TO TOT
REPLACE ALL REEND WITH 1
MARK1 = 1
SW2=0
DO WHILE SW2=0 ROLL
  IF MARK1 >= TOT
  PACK
  COUNT TO UNIQUE
  COUNT TO NEWGENES FOR D='H'.OR.D='O'
  5₩2=1
  LOOP
  ENDIF
GO MARKI
DUP - 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTE
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE REEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDEO ROLL
GO TOP
STORE R TO FUNC
USE "Analysis function.dbf"
LOCATE FOR P=FUNC
REFLACE CLONES WITH TOT
REPLACE GENES WITH UNIQUE
REPLACE NEW WITH NEWGENES.
USE TEMP1
SORT ON RFEND/D TO TEMP2
USE TEMP2
SET HEADING ON
?? STR(UNICUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
?? ' clones'
                         V Coincidence'
list off fields number, RFEND, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR, LENGTH, INIT, I
*SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 FIXELS FORT "Geneva",12 COLOR 0,0,
*list off fields RFEND, S, DESCRIPTOR
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE TEMPDESIG
```

PCT/US95/01160

```
* CCMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
USB TEMP1
COUNT TO TOT
REPLACE ALL REEND WITH 1
MARK1 = 1
SW2=0
DO WHILE SW2=0 ROLL
  IF MARK1 >= TOT
  PACK
  COUNT TO UNIQUE
  SW2=1
  LOOP
  ENDIF
GO MARKI
DUP = 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTS
  IF TESTA - TESTS
  DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
5W=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
GO TOP
STORE F TO DIST
USE "Analysis distribution.dbf"
LOCATE FOR P-DIST
REPLACE CLONES WITH TOT
REPLACE GENES WITH UNIQUE
USE TEMP1
sort on rfend/d to TEMP2
USE TEMP2
?? STR(UNIQUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
?? 'clones'
? '
                            V Coincidence'
list off fields number, RFEND, L. D. F. Z. R. C. ENTRY, S. DESCRIPTOR, LEWGIH, INIT, I
*SET PRINT OFF
CLOSE DATABASES
BRASE TEMP1.DBF
ERASE TEMP2, DBF
USB TEMPDESIG
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 USE TEMP1
 COUNT TO TOT
REPLACE ALL REEND WITH 1
 MARK1 = 1
 SW2=0
 DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
   SW2=1
   LOOP
  ENDIF
 GO MARK1
DUP = 1
 STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
 SKIP
 STORE ENTRY TO TESTE
  IF TESTA = TESTB
DELETE
  DUP .= DUP+1
  LOOP
  ENDIF
GO MARK1
REPLACE - RFEND WITH DUP
MARK1 = MARK1+DUP
SW=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
GO TOP
USE TEMP1
?? STR(UNIQUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
?? ' clones'
                            V Coincidence'
list off fields number, RFEND, L, D, F, Z, R, C, EMTRY, S, DESCRIPTOR, LENGTH, INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DRF
USE TEMPDESIG
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 USE "SmartGuy: FoxBASE+/Mac: fox files: Clomes.dbf"
 COPY TO TEMP1 FOR
 USE TEMP1
 COUNT TO IDGENE FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='R'.OR.D='A'
 DELETE FOR D='N'.OR.D='D'.OR.D='A'.OR.D='U'.OR.D='S'.OR.D='M'.OR.D='R'.OR.D='V'
 PACK
 COUNT TO TOT
 REPLACE ALL REEND WITH 1
 MARK1 = 1
 SW2=0
 DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
   SW2=1
   LOOP
   ENDIF
 GO MARKI
DUP = 1
STORE ENTRY TO TESTA
5W = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTE
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
  LOOP
  ENDIF
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
6W=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON RFEND/D, NUMBER TO TEMP2
USE TEMP2
REPLACE ALL START WITH RFEND/IDGENE*10000
?? STR(UNIQUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
?? ' clones'
     Coincidence V
                          V Clomes/10000'
set heading off
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
list fields number, RFEND, START, L, D, F, Z, R, C, EVTRY, S, DESCRIPTOR, INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1. DBP
ERASE TEMP2. DBF
USE *SmartGry:FoxBASE+/Mac:fox files:clones.dbf*
```

```
* COMPRESSION SUBROUTINE FOR ANALYSIS PROGRAMS
 COUNT TO IDGENE FOR D='E'.OR.D='O'.OR.D='H'.OR.D='N'.OR.D='R'.OR.D='A'
 DELETE FOR D='N'.OR.D='D'.OR.D='A'.OR.D='U',OR.D='S'.OR.D='M'.OR.D='R'.OR.D='V'
 COUNT TO TOT
 REPLACE ALL RFEND WITH 1
 MARK1 = 1
 SW2=0
 DO WHILE SW2=0 ROLL
   IF MARK1 >= TOT
   PACK
   COUNT TO UNIQUE
  SW2=1
  LOOP
  ENDIF
GO MARKI
DUP = 1
STORE ENTRY TO TESTA
SW = 0
DO WHILE SW=0 TEST
SKIP
STORE ENTRY TO TESTE
  IF TESTA = TESTB
  DELETE
  DUP = DUP+1
LOOP
  ENDIF
GO MARKI
REPLACE RFEND WITH DUP
MARK1 = MARK1+DUP
5W=1
LOOP
ENDDO TEST
LOOP
ENDDO ROLL
*BROWSE
*SET PRINTER ON
SORT ON RFEND/D, NUMBER TO TEMP2
USE TEMP2
REPLACE ALL START WITH RFEND/IDGENE*10000
?? STR(UNIQUE,5,0)
?? 'genes, for a total of '
?? STR(TOT,5,0)
?? ' clones'
? ' Coincidence V
                           V Clones/10000'
set heading off
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0,
list fields number, RFEND, START, L, D, F, Z, R, C, ENTRY, S, DESCRIPTOR; INIT, I
*SET PRINT OFF
CLOSE DATABASES
ERASE TEMP1.DBF
ERASE TEMP2.DBF
USE *SmartGuy: FoxBASE+/Mac: fox files:clones.dbf*
```

```
USE TEMP1
COUNT TO TOT
?? 'Total of'
?? STR(TOT,4,0)
?? 'clones'
?
*list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR,LENGTH,RFEND,INIT,I
list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR
CLOSE DATABASES
ERASE TEMP1.DBF
USE TEMPDESIG
```

```
*Lifescan memu; version 8-7-94
SET TALK OFF
set device to screen
CLEAR
USE "SmartGuy: FoxBASE+/Mac:fox files:clones.dbf"
STORE LUPDATE() TO Update
GO BOTTOM
STORE RECNO() TO cloneno
STORE 6 TO Chooser
DO WHILE .T.
* Program.: Lifeseg menu.fmt
  Date...: 1/11/95
  Version .: FootBASE+/Mac, revision 1.10
* Notes....: Format file Lifeseg menu
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",268 COLOR 0,0,
@ PIXELS 18,126 TO 77,365 STYLE 28479 COLOR 32767,-25600,-1,-16223,-16721,-15725 B PIXELS 110,29 TO 188,217 STYLE 3871 COLOR 0,0,-1,-25600,-1,-1
PIXELS 45,161 SAY "LIFESEQ" STYLE 65536 FONT "Geneva",536 COLOR 0.0,-1,-1,7135,5884
PIXELS 36,269 SAY "TM" STYLE 65536 FONT "Geneva",12 COLOR 0.0,-1,-1,7135,5884
PIXELS 63,143 SAY "Molecular Biology Desktop" STYLE 65536 FONT "Helvetica",18 COLOR 0.0.0,
@ PIXELS 90,252 TO 251,467 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1
@ PIXELS 45,296 SAY "V1.30" STYLE 65536 FONT "Geneva",782 COLOR 0,0,-1,-1,-1,-1
* EOF: Lifeseq menu.fmt
READ
DO CASE
CASE Chooser=1
DO "SmartGuy:FoxPASE+/Mac:fox files:Output programs:Naster analysis 3.prg"
CASE Chooser=2
DO "SmartGuy: Fox BASE+/Mac: fox files: Output programs: Subtraction 2.prg"
CASE Chooser=3
DO "SmartGuy:FoxBASE+/Mac:fox files:Output programs:Northern (single).prg"
CASE Chooser=4
USE "Libraries.dbf"
BROWSB
CASE Chooser=5
DO "SmartGuy: FoxEASE+/Mac: fox files: Output programs: See individual clone.prg"
DO "SmartGuy:FoxBASE+/Mac:fox files:Libraries:Output programs:Menu.prg"
CASE Chooser=7
CLEAR
SCREEN 1 OFF
RETURN
ENDCASE
LOOP
```

ENDDO

PCT/US95/01160

? 'Arranged by FUNCTION'

```
61,30 SAY "Database Subset Analysis" STYLE 65536 FONT "Geneva",274 COLOR 0,0,0,-1,-1,-1
? date()
?? TIMB()
? 'Clone numbers '
?? STR (INITIATE, 6,0)
?? 'through '
?? STR (TERMINATE, 6, 0)
? Libraries:
IP ENTIRE=1
? 'All libraries'
ENDIF
IF ENTIRE=2
    MARK=1
    DO WHILE .T.
   IF MARK>STOPIT
    EXIT
    ENDIF
    USE SELECTED
    GO MARK
    ?? TRIM(libname)
    STORE MARK+1 TO MARK
    LOOP
    ENDDO
ENDIF
? 'Designations: '
IF Ematch=0 .AND, Hmatch=0 .AND, Omatch=0
?? 'All'
ENDIF
IF Ematch=1
?? 'Exact,'
ENDIF
IF Hmatch=1
?? 'Human, '
ENDIF
IF Omatch=1
?? 'Other sp.'
ENDIF
IF CONDEN=1
? 'Condensed format analysis'
ENDIF
IF ANAL=1
? 'Sorted by NUMBER'
ENDIF
IF ANAL=2
? 'Sorted by ENTRY'
ENDIF
IF ANAL-3
? 'Arranged by ABUNDANCE'
ENDIF
IF ANAL=4
? 'Sorted by INTEREST'
ENDIF
IP ANAL=5
? 'Arranged by LOCATION'
ENDIF
IF ANAL-5
? 'Arranged by DISTRIBUTION'
ENDIF
IF ANAL=7
```

```
PNDIF
? 'Total clones represented: '
?? STR(STARTOT,6.0)
? 'Total clones analyzed: '
?? STR(ANALTOT,6.0)
?
```

```
USE TEMP1
COUNT TO TOT
?? 'Total of'
?? STR(TOT,4,0)
?? 'clones'
?
*list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR,LENGTH,RFEND,INIT,I
list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR
CLOSE DATABASES
ERASE TEMP1.DBF
USE TEMPDESIG
```

```
USE TEMP1
COUNT TO TOT
?? ' Total of'
?? STR(TOT,4,0)
?? ' clones!
?
*list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR,LENGTH,RFEND,INIT,I
list off fields number,L,D,F,Z,R,C,ENTRY,DESCRIPTOR
CLOSE DATABASES
ERASE TEMP1.DB?
USE TEMPDESIG
```

```
*Northern (single), version 11-25-94
 close databases
 SET TALK OFF
 SET PRINT OFF
 SET EXACT OFF
 CLEAR
 STORE
                                ' TO Eobject
 STORE '
                                                                    ' TO Dobject
 STORE 0 TO Numb
 STORE 0 TO Zog
STORE 1 TO Bail
DO WHILE .T.
 * Program.: Northern (single).fmt
* Date...: 8/8/94
 * Version .: FoxBASE+/Mag, revision 1.10
 * Notes...: Format file Northern (single)
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",12 COLOR 0,0,0 @ PIXELS 15,81 TO 46,397 STYLE 28447 COLOR 0,0,-1,-25600,-1,-1 @ PIXELS 89,79 TO 192,422 STYLE 28447 COLOR 0,0,0,-25600,-1,-1 @ PIXELS 115,98 SAY "Entry #: " STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
PIXELS 115,98 SAY "Entry style 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
@ PIXELS 115.173 GET Eobject STYLE 0 FONT "Geneva",12 SIZE 15.142 COLOR 0,0,0,-1,-1,-1
@ PIXELS 145.89 SAY "Description" STYLE 65536 FONT "Geneva",12 COLOR 0,0,0,-1,-1,-1
@ PIXELS 145,173 GET Dobject STYLE 0 FONT "Geneva",12 SIZE 15,241 COLOR 0,0,0,-1,-1,-1
@ PIXELS 35,89 SAY "Single Northern search screen" STYLE 65536 FONT "Geneva",274 COLOR 0,0,-
@ PIXELS 220,162 GET Bail STYLE 65536 FONT "Chicago",12 PICTURE "@*R Continue;Bail out" SIZE
@ PIXELS 175,98 SAY "Clone #:" STYLE 65536 FONT "Geneva";12 COLOR 0,0,0,-1,-1,-1
@ PIXELS 175,173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,0,-1,-1,-1
@ PIXELS 175,173 GET Numb STYLE 0 FONT "Geneva",12 SIZE 15,70 COLOR 0,0,0,-1,-1,-1
@ PIXELS 80,152 SAY "Enter any ONE of the following: STYLE 65536 FONT "Geneva",12 COLOR -1,
* EOF: Northern (single).fmt
READ
IF Bail=2
CLEAR
screen 1 off
RETURN
EMDIF
USE "SmartGuy: FoxBASE+/Mac:Fox files:Lookup.dbf"
SET TALK ON
IF Eobject<>'
STORE UPPER (Eobject) to Eobject
SET SAFETY OFF
SORT ON Entry TO "Lookup entry.dbf"
SET SAFETY ON
USE "Lookup entry.dbf"
LOCATE FOR Look=Eobject
IF .NOT.FOUND()
CLEAR
LOOP
ENDIF
BROWSE
STORE Entry TO Searchval
CLOSE DATABASES
ERASE "Lookup entry.dbf"
ENDIP
IF Dobject<>'
SET EXACT OFF
SET SAFETY OFF
SORT ON descriptor TO "Lookup descriptor.dbf"
SET SAFETY On
USB "Lookup descriptor.dbf"
LOCATE FOR UPPER (TRIM (descriptor)) = UPPER (TRIM (Dobject))
IF .NOT.FOUND()
CLEAR
```

```
LOOP
 ENDIF
BROWSE
STORE Entry TO Searchval
CLOSE DATABASES
ERASE "Lookup descriptor.dbf"
 SET EXACT ON
ENDIP
IF Numb<>0
USE "SmartGuy: FoxBASE+/Mac: Fox files: clones.dbf"
GO Numb
BROWSE
STORE Entry TO Searchval
ENDIP
CLEAR
? 'Northern analysis for entry '
?? Searchval
? 'Enter Y to proceed'
NO OT TIAW
CLEAR
IF UPPER (OK) <> 'Y'
screen 1 off
RETURN
ENDIF
* COMPRESSION SUBROUTINE FOR Library, dbf
7 'Compressing the Libraries file now...'
USE 'SmartGuy:FoxBASE+/Mac:Fox files:libraries.dbf'
SET SAFETY OFF
SORT ON library TO "Compressed libraries.dbf"
* FOR entered>0
SET SAFETY ON
USE "Compressed libraries.dbf"
DELETE FOR entered=0
PACK
COUNT TO TOT
MARK1 = 1
$₩2±0
DO WHILE SW2=0 ROLL
  IF MARK1 >= TOT
  PACK
  5W2=1
  LOOP
  ENDIP
GO MARKI
STORE library TO TESTA
SKIP
STORE Library TO TESTS
IF TESTA = TESTS
DELETE
ENDIF
MARK1 = MARK1+1
LOOP
ENDDO ROLL
* Northern analysis
CLEAR
? 'Doing the northern now...'
USE "SmartGuy: FoxBASE+/Mac: Fox files: clones.dbf"
SET SAFETY OFF
COPY TO "Hit .dbf" FOR entry=searchval
SET SAFETY ON
```

PCT/US95/01160

```
CLOSE DATABASES
SELECT 1
USE "Compressed libraries.dbf"
STORE RECCOUNT() TO Entries
SELECT 2
USE "Hits.dbf"
Mark=1
DO WHILE .T.
SELECT 1
IF MarkoEntries
EXIT
ENDIF
GO MARK
STORE library TO Jigger
SELECT 2
COUNT TO Zog FOR library=Jigger
SELECT 1
REPLACE hits with Zog
Mark=Mark+1
LOOP
ENIDDO .
SELECT 1
BROWSE FIELDS LIMRARY, LIENAME, ENTERED, HITS AT 0,0
CLEAR
? 'Enter Y to print:'
IF UPPER (PRINSET) = 'Y'
SET PRINT ON
CLEAR
PJECT.
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FONT "Geneva",14 COLOR 0,0,0 ? 'DATABASE ENTRIES MATCHING ENTRY '
?? Searchval
? DATE()
SCREEN 1 TYPE 0 HEADING "Screen 1" AT 40,2 SIZE 286,492 PIXELS FORT "Geneva",7 COLOR 0,0,0, LIST OFF FIELDS library, librare, entered, bits
SELECT 2
LIST OFF FIELDS NUMBER, LIBRARY, D.S.F.Z.R. ENTRY, DESCRIPTOR, RFSTART, START, RFEND
SET TALK OFF
SET PRINT OFF
ENDIP
CLOSE DATABASES
SET TALK OFF
CLEAR
DO "Test print.prg"
RETURN
```

· TABLE 6

libname library Inflamed adenoid ADENINB01 ADRENORO1 Adrenal gland (1)
ADRENOTO1 Adrenal gland (T) AMLBNOTO1 AML blast cells (1) BMARNOTO1 Bone marrow BMARNOTO2 Bone marrow (T)
CARDNOTO1 Cardiac muscle (T) CHAONOTO1 Chin. hamster ovary CORNNOTO1 Comeal stroma FIBRAGT01 Fibroblast, AT 5 FIBRAGT02 Fibroblast, AT 30 FIBRANTO1 Fibroblast, AT FIERNGT01 Fibroblast, uv 5 FIBRNGT02 Fibroblast, uv 30 FIERNOTO: Fibroblast FIERNOTO2 Fibroblast, normal HMC1NOT01 Mast cell line HMC-1 HUVELPBO1 HUVECIFN,TNF,LPS HUVENOBO1 HUVEC control **HUVESTBO1 HUVEC** shear stress HYPONKIDNNOTO1
LIVRNOTO1
LUNGNOTO1
MUSCNOTO1
OVIDNOB01
PANCNOTO1
PITUNOR01
PITUNOR01
PITUNOR01
PItuliary (r)
Placenta
Placenta
Placenta HYPONOB01 Hypothelamus SINTNOTO2 SPLNFETO1 SPLNNOTO2 Small intestine (T) Spieen+liver, fetal Spieen (T) Stornach Rheum, synovium STOMNOTO SYNORAB01 TBLYNOT01 T + B lymphoblast TESTNOTO1 Testie (T) THP1NOB01 THP-1 control THP1PEB01 THP phorbol THP-1 phorbol LPS THP1PLB01 U937NOT01 U937, monocytic leuk

numberlibrary	dsizte	entry	descriptor	rfstari	81811	rfend
2304 U837NOTO	SEHCCT	HUMEF18	Elongation lactor 1-beta	₽-	0	773
3240 HMC1NOTO	1 EHCCT	HUMEFIB	Elongation (actor 1-beta	٥	370	773
3259 HMC1NOTO	LEHCCT	HUMEF18	Elongation factor 1-beta	0	371	773
	1 EHCCT		Elongation factor 1-beta	0	470	773
	EHCCT		Elongation factor 1-beta	٥	327	773
	EHCCT		Elongation factor 1-beta	. 0	375	773

WHAT IS CLAIMED IS:

5

1. A method of analyzing a specimen containing gen transcripts, said method comprising the steps of:

- (a) producing a library of biological sequences;
- (b) generating a set of transcript sequences, where each of the transcript sequences in said set is indicative of a different one of the biological sequences of the library;
- (c) processing the transcript sequences in a

 10 programmed computer in which a database of reference
 transcript sequences indicative of reference biological
 sequences is stored, to generate an identified sequence
 value for each of the transcript sequences, where each said
 identified sequence value is indicative of a sequence

 15 annotation and a degree of match between one of the
 transcript sequences and at least one of the reference
 transcript sequences; and
 - (d) processing each said identified sequence value to generate final data values indicative of a number of times0 each identified sequence value is present in the library.
 - 2. The method of claim 1, wherein step (a) includes the steps of:

obtaining a mixture of mRNA; making cDNA copies of the mRNA;

- isolating a representative population of clones transfected with the cDNA and producing therefrom the library of biological sequences.
 - 3. The method of claim 1, wherein the biological sequences are cDNA sequences.
- 4. The method of claim 1, wherein the biological sequences are RNA sequences.
 - 5. The method of claim 1, wherein the biological sequences are protein sequences.

6. The method of claim 1, wherein a first value of said degree of match is indicative of an exact match, and a second value of said degree of match is indicative of a non-exact match.

- 7. A method of comparing two specimens containing gene transcripts, said method comprising:
 - (a) analyzing a first specimen according to the method of claim 1;
- (b) producing a second library of biological10 sequences;
 - (c) generating a second set of transcript sequences, where each of the transcript sequences in said second set is indicative of a different one of the biological sequences of the second library;
- (d) processing the second set of transcript sequences in said programmed computer to generate a second set of identified sequence values known as further identified sequence values, where each of the further identified sequence values is indicative of a sequence annotation and a degree of match between one of the biological sequences of the second library and at least one of the reference sequences;
- (e) processing each said further identified sequence value to generate further final data values indicative of a
 25 number of times each further identified sequence value is present in the second library; and
- (f) processing the final data values from the first specimen and the further identified sequence values from the second specimen to generate ratios of transcript
 sequences, each of said ratio values indicative of differences in numbers of gene transcripts between the two specimens.
- 8. A method of quantifying relative abundance of mRNA in a biological specimen, said method comprising the steps35 of:
 - (a) isolating a population of mRNA transcripts from the biological specimen;

(b) identifying g nes from which the mRNA was transcribed by a sequence-specific method;

- (c) determining numbers of mRNA transcripts c rresponding to each of the genes; and
- 5 (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts.
 - 9. A diagnostic method which comprises producing a gene transcript image, said method comprising the steps of:
- 10 (a) isolating a population of mRNA transcripts from a biological specimen;
 - (b) identifying genes from which the mRNA was transcribed by a sequence-specific method;
- (c) determining numbers of mRNA transcripts
 15 corresponding to each of the genes; and
- (d) using the mRNA transcript numbers to determine the relative abundance of mRNA transcripts within the population of mRNA transcripts, where data determining the relative abundance values of mRNA transcripts is the gene transcript image of the biological specimen.
 - 10. The method of claim 9, further comprising:
 - (e) providing a set of standard normal and diseased gene transcript images; and
- (f) comparing the gene transcript image of the
 25 biological specimen with the gene transcript images of step
 (e) to identify at least one of the standard gene
 transcript images which most closely approximate the gene
 transcript image of the biological specimen.
- 11. The method of claim 9, wherein the biological 30 specimen is biopsy tissue, sputum, blood or urine.
 - 12. A method of producing a gene transcript image, said method comprising the steps of
 - (a) obtaining a mixture of mRNA;
 - (b) making cDNA copies of the mRNA;

(c) inserting the cDNA into a suitable vector and using said vector to transfect suitable host strain cells which are plated out and permitted to grow into clones, each clone representing a unique mRNA;

- (d) isolating a representative population of recombinant clones;
 - (e) identifying amplified cDNAs from each clone in the population by a sequence-specific method which identifies gene from which the unique mRNA was transcribed;
- 10 (f) determining a number of times each gene is represented within the population of clones as an indication of relative abundance; and
- (g) listing the genes and their relative abundance in order of abundance, thereby producing the gene transcript15 image.
 - 13. The method of claim 12, also including the step of diagnosing disease by:

repeating steps (a) through (g) on biological specimens from random sample of normal and diseased humans, encompassing a variety of diseases, to produce reference sets of normal and diseased gene transcript images;

obtaining a test specimen from a human, and producing a test gene transcript image by performing steps (a) through (g) on said test specimen;

comparing the test gene transcript image with the reference sets of gene transcript images; and

identifying at least one of the reference gene transcript images which most closely approximates the test gene transcript image.

30 14. A computer system for analyzing a library of biological sequences, said system including:

35

means for receiving a set of transcript sequences, where each of the transcript sequences is indicative of a different one of the biological sequences of the library; and

means for processing the transcript sequences in the computer system in which a database of reference transcript

sequences indicative of reference biological sequences is stored, wherein the comput r is programmed with software for generating an identified sequence value for each of the transcript sequences, where each said identified sequence value is indicative of a sequence annotation and a degree of match between a different one of the biological sequences of the library and at least one of the reference transcript sequences, and for processing each said identified sequence value to generate final data values indicative of a number of times each identified sequence value is present in the library.

- 15. The system of claim 14, also including:
 library generation means for producing the library of
 biological sequences and generating said set of transcript
 sequences from said library.
 - 16. The system of claim 15, wherein the library generation means includes:

means for obtaining a mixture of mRNA; means for making cDNA copies of the mRNA;

20

means for inserting the cDNA copies into cells and permitting the cells to grow into clones;

means for isolating a representative population of the clones and producing therefrom the library of biological sequences.

SYBASE database Structure Library Preparation

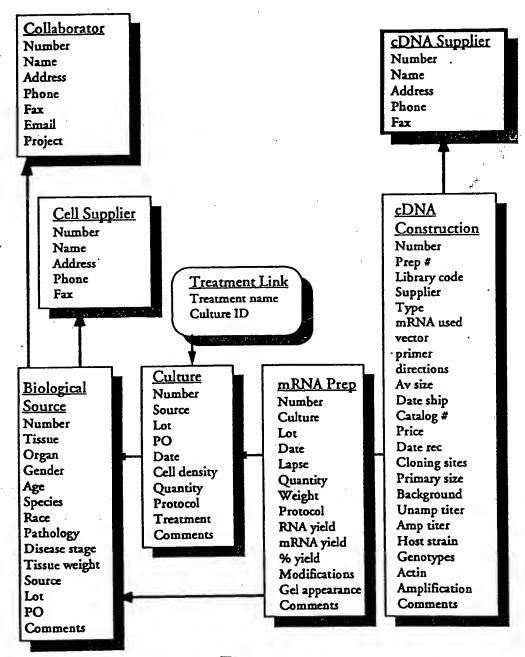


Figure 1

Identified Sequences Tempnum Templesig Tempsub Templib Libsort Trosduz Tempsubsost Temptorsort Final Data

Figure 2

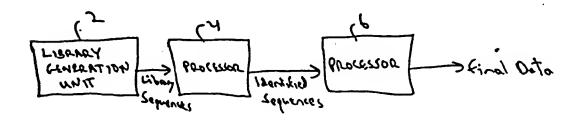


Figure 3

Incyte Bioinformatics Process

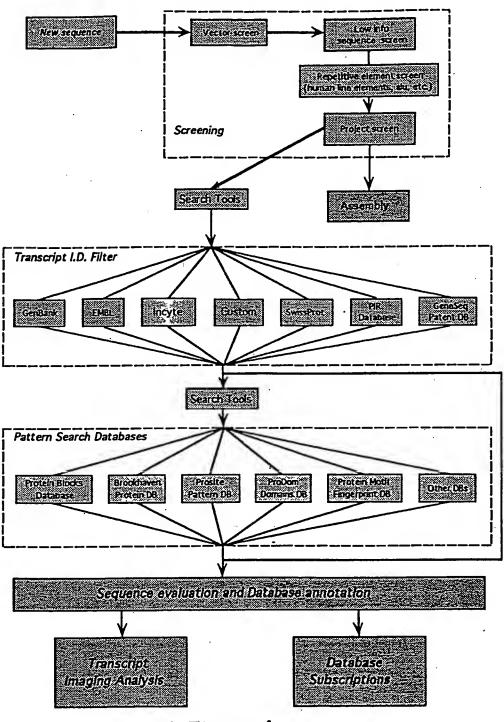


Figure 4

INTERNATIONAL SEARCH REPORT

International application No. PCT/US95/01160

A. CLASSIFICATION OF SUBJECT MATTER							
IPC(6) :C12Q 1/68; G06F 15/00							
US CL: 435/6; 364/413.02 According to International Patent Classification (IPC) or to both national classification and IPC							
	LDS SEARCHED						
	ocumentation searched (classification system followed	by classification symbols)					
U.S. :	435/6; 364/413.02	•					
Documental	tion searched other than minimum documentation to the	extent that such documents are included	in the fields searched				
			' '				
Electronic o	lata base consulted during the international search (na	ne of data base and, where practicable,	, search terms used)				
	UNE, APS, transcript, transcripts, cdan#, mmai						
CAS ON	Disc, Al O, Unitedity Childson Du, Courte, Tillion						
C. DOC	CUMENTS CONSIDERED TO BE RELEVANT						
Category*	Citation of document, with indication, where ap	propriate, of the relevant passages	Relevant to claim No.				
X	IntelliGenetics Suite, Release 5.4, A		15 and 16				
	issued January 1993 by IntelliGe						
Υ	Camino Real, Mountain View, C		1-14				
	States of America, pages (1-6)-(1-	19) and (2-9)-(2-14), see	·				
	entire document.						
			4.40				
Υ	Science, Volume 252, issued 21 Ju		1-16				
	al, "Complementary DNA sequenc						
	tags and human genome project"	, pages 1651-1656, see					
	entire document.						
		<u>.</u>					
	·						
-							
X Further documents are listed in the continuation of Box C. See patent family annex.							
Special categories of cited documents: T later document published after the international filing date or priority date and not in conflict with the application but cited to understand the							
"A" document defining the general state of the art which is not considered principle or theory underlying the invention to be of particular relevance							
"E" carrier document published on or after the international filing date "X" document of particular relevance; the claimed invention cannot be considered povel or cannot be considered to involve an inventive step							
L document which may throw doubte on priority claim(s) or which is when the document is taken alone							
cited to establish the publication date of smother citation or other special reason (as specified) document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is							
"O" document referring to an oral disclosure, use, exhibition or other combined with one or more other such documents, such combination being obvious to a person skilled in the art							
"P" document published prior to the international filing date but later than "&" document member of the same patent family							
Date of the actual completion of the international search Date of mailing of the international search report							
0.4 MAY 1995							
27 APRII	L 1995	0 I III					
Name and mailing address of the ISA/US Authorized officer							
Name and mailing address of the ISA/US Commissioner of Patents and Trademarks Box PCT Westington D.C. 20031							
Washingto	wantiguit p.c. pap.						
English: 1- 1	v. (202) 205_2220	Telephone No. (703) 308-0196					

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US95/01160

C (Continua	tion). DOCUMENTS CONSIDERED TO BE RELEVANT			
Category*	Citation of document, with indication, where appropriate, of the releva	ant passages	Relevant to claim No.	
Y	Nucleic Acids Research, Volume 19, No. 25, issued 19 Hara et al, "Subtractive cDNA cloning using oligo(dT) ₃ PCR: isolation of cDNA clones specific to undifferenti embryonal carcinoma cells", pages 7097-7104, see entire document.	1-16		
X Y	Nature Genetics, Volume 2, No. 3, issued November 1 Okubo et al, "Large scale cDNA sequencing for analysi quantitative and qualitative aspects of gene expression", 173-179, see narrative text portion of entire document.	1, 3 2 and 4-16		
ì				
`,				
	·			
	·			
	,			
·	·			
			·	